



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
STATE PUBLIC HEALTH LABORATORY
BREATH ALCOHOL PROGRAM
INTOX EC/IR II MAINTENANCE REPORT

REPORT #3

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days). Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service. Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

INTOX EC/IR II SN 12814	NAME OF AGENCY MARYLAND HEIGHTS POLICE	DATE OF INSPECTION 08/26/2016
----------------------------	---	----------------------------------

LOCATION OF INSTRUMENT (STREET AND CITY) 11911 DORSETT RD MARYLAND HEIGHTS	TIME OF INSPECTION 05:50 CDT
---	---------------------------------

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or is operating within established limits. (Write in observed values where determined). Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC RECORD	
<input checked="" type="checkbox"/> BLANK CHECK	<input checked="" type="checkbox"/> CO2 CHECK
<input checked="" type="checkbox"/> FC 1 TEMP	<input checked="" type="checkbox"/> FLOW CHECK
<input checked="" type="checkbox"/> SRC TEMP	<input checked="" type="checkbox"/> FCB CHECK
<input checked="" type="checkbox"/> DET TEMP	<input checked="" type="checkbox"/> CRC COMP CHECK
<input checked="" type="checkbox"/> BT TEMP	<input checked="" type="checkbox"/> CRC CAL CHECK
<input checked="" type="checkbox"/> STD 2 TEMP	<input checked="" type="checkbox"/> PRINT TEST
<input checked="" type="checkbox"/> ETH CHECK	

BREATH ANALYZER ACCURACY STANDARDS	
<input checked="" type="checkbox"/> SIMULATOR SOLUTION	<input type="checkbox"/> COMPRESSED ETHANOL-GAS MIXTURE

<input checked="" type="checkbox"/> STANDARD SUPPLIER	GUTH LABORATORIES	LOT# 15050	EXP. DATE 03/09/2017
---	-------------------	------------	----------------------

<input checked="" type="checkbox"/> SIMULATOR TEMP (34°C ±0.2°C) 34°C +/- .2°	SIMULATOR S/N SD2240	SIMULATOR EXP DATE 10/13/2016
--	-------------------------	----------------------------------

CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within +5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

0.10% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE

0.08% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE

0.04% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 0.096 g/210L	TEST 2 0.097 g/210L	TEST 3 0.097 g/210L
---------------------	---------------------	---------------------

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:

REFUSALS	0	0-.04	36	.05-.09	0	.10-.14	0	.15-.19	0	OVER .19	0
----------	---	-------	----	---------	---	---------	---	---------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

UNIT & SIMULATOR MEET DHSS SPECIFICATIONS. DMS230

INSPECTING OFFICER	
SIGNATURE 	PRINT FULL NAME STOLARSKI
TYPE II PERMIT NUMBER 250056	EXPIRATION DATE 02/20/2017
	TELEPHONE NUMBER (314) 298-8700

RETURN COMPLETED REPORT TO THE:
Breath Alcohol Program, Missouri Department of Health and Senior Services,
Southeast District Office, 2875 James Blvd, Poplar Bluff, MO 63901



GUTH LABORATORIES, INC.

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-664-6470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **15050** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **March 11, 2015**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.1218%** (w/vol) ethyl alcohol. The expiration date for this lot number is **March 9, 2017** at 11:59 PM.

When used in a calibrated Simulator, operating at $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$, this solution will give a breath alcohol analysis instrument reading of **0.100 g/210L \pm 3%**.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN08051301 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.



STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

2

PERMIT
TYPE II

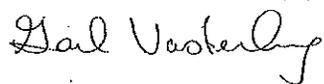
DANIEL M STOLARSKI

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

INTOX EC/IR II

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 2/20/2015
 NUMBER 250056
 EXPIRES 2/20/2017


 DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

 ,acting director
 DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

MO 580-0771 (6-10)

LAB-4 (R6-10)

 STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

INSTRUMENT OPERATOR CARD

The named cardholder is authorized to operate an evidential breath alcohol instrument for the determination of the alcoholic content in breath form of expired air in Missouri.



Operator STOLARSKI, DANIEL
 Permit No 250056
 Date issued 2/20/2015 Date Expires 2/20/2017