



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
 STATE PUBLIC HEALTH LABORATORY
ALCO-SENSOR IV WITH PRINTER MAINTENANCE REPORT

RECEIVED
 By Carol Day at 11:40 am, Jul 14, 2014

Complete this report in duplicate at the time of the regular monthly preventative maintenance check, and whenever instrument is repaired. Send copy to Department of Health and Senior Services; retain original in department file.

ALCO SENSOR IV SN 097403	PRINTER SN 096.3580.877	DATE OF INSPECTION 07/10/2014
LOCATION OF INSTRUMENT (STREET AND CITY) 317 BROADWAY, ELSBERRY		TIME OF INSPECTION 12:20 pm

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

DIGITAL READOUT (ALL ELEMENTS OPERATIONAL)

TEMPERATURE OF ALCO SENSOR (10°C - 40°C)

PRINTER WORKING PROPERLY

TIME AND DATE DISPLAYING PROPERLY

BREATH ALCOHOL ACCURACY STANDARDS

SIMULATOR SOLUTION COMPRESSED ETHANOL-GAS MIXTURE

STANDARD SUPPLIER GUTH LABS LOT # 14110 EXP. DATE 05/01/2016

SIMULATOR TEMPERATURE (34°C ± 0.2°C) 34.0 SIMULATOR SN SD2287 SIMULATOR EXP DATE 02/10/2015

CALIBRATION CHECK – (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Check the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

0.100% STANDARD - MUST READ BETWEEN 0.095% and 0.105% INCLUSIVE

0.080% STANDARD - MUST READ BETWEEN 0.076% and 0.084% INCLUSIVE

0.040% STANDARD - MUST READ BETWEEN 0.038% and 0.042% INCLUSIVE

TEST 1 ➔ .101

TEST 2 ➔ .102

TEST 3 ➔ .102

RFI DETECTOR OPERATING

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT: (DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	1	(0-.04)	0	(.05-.09)	0	(.10-.14)	0	(.15-.19)	0	(OVER .19)	0
----------	---	---------	---	-----------	---	-----------	---	-----------	---	------------	---

List any new parts and describe any alteration or modification that was made to restore the instrument to operate satisfactorily and within established limits (use other side if necessary).

INSPECTING OFFICER

SIGNATURE

PRINT NAME
 ROBERT J BODLEY

TYPE II PERMIT NUMBER/EXPIRATION DATE
 240121 03/24/2016

TELEPHONE NUMBER
 (573) 898-5456

Return completed report to the: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
 2875 James Boulevard
 Poplar Bluff, MO 63901

AS IV Serial no: 097403
Version no: 004C

TEST RECORD 00679

Temp Date Time ^{s/} 210L

Air Blank:
07/10/14 12:20 .000
Calibration Check:
24 07/10/14 12:20 .101

Subject Name

Subject I.D.

Operator Name, I.D.

Location

AS IV Serial no: 097403
Version no: 004C

TEST RECORD 00681

Temp Date Time ^{s/} 210L

Air Blank:
07/10/14 12:23 .000
Calibration Check:
24 07/10/14 12:23 .102

Subject Name

Subject I.D.

Operator Name, I.D.

Location

AS IV Serial no: 097403
Version no: 004C

TEST RECORD 00680

Temp Date Time ^{s/} 210L

Air Blank:
07/10/14 12:21 .000
Calibration Check:
24 07/10/14 12:21 .102

Subject Name

Subject I.D.

Operator Name, I.D.

Location

AS IV Serial no: 097403
Version no: 004C

TEST RECORD 00682

Temp Date Time ^{s/} 210L

Void: RFI
12 07/10/14 12:24

Subject Name

Subject I.D.

Operator Name, I.D.

Location



GUTH LABORATORIES, INC.

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **14110** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **May 5, 2014**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.1206%** (w/vol) ethyl alcohol. The expiration date for this lot number is **May 1, 2016** at 11:59 PM.

When used in a calibrated Simulator, operating at $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$, this solution will give a breath alcohol analysis instrument reading of **0.100 g/210L \pm 3%**.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN122211-02 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.



STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

2

**PERMIT
 TYPE II**

ROBERT J BODLEY

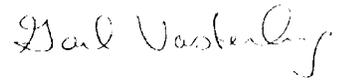
is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

ALCO-SENSOR IV WITH PRINTER

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE **3/24/2014**
 NUMBER **240121**
 EXPIRES **3/24/2016**


 DIRECTOR OF STATE PUBLIC HEALTH LABORATORY


 DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES