



STATE PUBLIC HEALTH LABORATORY
BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

RECEIVED
By Carol Day at 11:29 am, Sep 08, 2014

REPORT #6

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).
Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.
Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 201182	NAME OF AGENCY Hollister Police Department	DATE OF INSPECTION 09/04/2014
LOCATION OF INSTRUMENT (STREET AND CITY) #1 Gage Drive, Hollister		TIME OF INSPECTION 0:10 am

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) <u>9/4/14, 0009</u>
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER <u>49</u> °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

INDICATOR LIGHTS

SIMULATOR SOLUTION SUPPLIER RepCo Marketing Inc. LOT # 14001 EXP. DATE 04/30/2016

SIMULATOR TEMP (34°C ± 0.2°C) 34.0 °C SIMULATOR SN SD2732 EXP. DATE 08/26/2015

CALIBRATION CHECK – (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE

0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE

0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 .101%	TEST 2 .100%	TEST 3 .100%
---------------	---------------	---------------

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:
(DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	3	(0-.04)	1	(.05-.09)	1	(.10-.14)	1	(.15-.19)	0	OVER .19	0
----------	---	---------	---	-----------	---	-----------	---	-----------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

Evidence Ticket Printer Unit Replaced and Tested, Works to Specifications

Instrument Operating within DOHSS Specs and guidelines

INSPECTING OFFICER	
SIGNATURE 	PRINT FULL NAME Sgt. Timothy E. Matthews
TYPE II PERMIT NUMBER/EXPIRATION DATE 220253 09/11/2014	TELEPHONE NUMBER (417) 334-6565

RETURN COMPLETED REPORT TO THE: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
2875 James Blvd.
Poplar Bluff, MO 63901

CERTIFICATE OF ANALYSIS

MANUFACTURER AND SUPPLIER: RepCo Marketing, Inc.
LOT NUMBER: 14001
EXPIRATION DATE: April 30, 2016 at 11:59 p.m.

RepCo Marketing, Inc. certifies the following:

RepCo Marketing, Inc. manufactured, tested and supplied Lot Number 14001 of Alcohol Certified Solution for simulators. Random samples of said lot number were analyzed by an independent laboratory utilizing a gas chromatograph and found to contain .1216 gms/dl +/- .003 gms/dl wt/vol ethanol (95% Confidence).

The alcohol and distilled water used in the solution were found to be free of any interfering substance.

This solution will produce a vapor alcohol value of .100 +/- 3% gms/210L Breath when heated to 34 Degrees Celsius +/- 0.2 Degrees Celsius in a simulator (95% Confidence).

The date of manufacture for this lot number is May 1, 2014
The expiration date for this lot number is April 30, 2016 at
11:59 p.m.

This document is a true representation of the original Certificate of Analysis.



Cecil B. Garner, President
RepCo Marketing, Inc.

Face This Side Down - This Edge In First

BAC DataMaster Evidence Ticket

STATE OF MISSOURI
HOLLISTER POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 201102
L 09/04/14

TESTING OFFICER:
MATTHEWS/T/E
OFFICER I.D.: 582
PERMIT NUMBER: 220252
EXPIRATION DATE: 09/11/14
MISCELLANEOUS DATA:

--- SUPERVISOR MODE ---

BLANK TEST	.000	00:13
INTERNAL STANDARD	VERIFIED	00:13
EXTERNAL STANDARD	.101	00:14
BLANK TEST	.000	00:15
EXTERNAL STANDARD	.100	00:15
BLANK TEST	.000	00:16
EXTERNAL STANDARD	.100	00:16
BLANK TEST	.000	00:17

N = 3
SIM. = .1
AVG. = .1003

Operator Signature



2208-02

Face This Side Down - This Edge In First

BAC DataMaster Evidence Ticket

STATE OF MISSOURI
HOLLISTER POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 201102
09/04/14
00109


--- DIAGNOSTIC CHECK ---

COMPUTER:	OKAY
PROGRAM (04-07-2009):	OKAY
HEATERS:	
SAMPLE CHAMBER:	495
FLOW DETECTOR:	OKAY
PUMP:	
HIGH SPEED:	OKAY
DETECTOR:	OKAY
FILTERS:	OKAY
QUARTZ STANDARD:	OKAY
CALIBRATION:	OKAY

PRINTER TEST

!@#\$%^&*+,-./0123456789:;<=>?@AB0CDEF
HIJKLMNOPQRSTUVWXYZ[\]^_`abcedfghijklmnop
qrstuvwxyz{|}~?>

Operator Signature

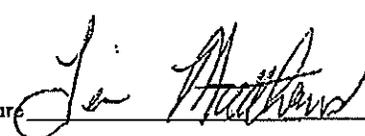

2208-02

Face This Side Down - This Edge In First

**BAC DataMaster
Evidence Ticket**

STATE OF MISSOURI
HOLLISTER POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 221120
09/04/14

ARREST TIME: 03:00
SUBJECT NAME:
TEST
DOB: 01/01/91 SEX: M
STATE/D.L.: MO/111
ARRESTING OFFICER:
TEST
OFFICER I.D.: 5
TESTING OFFICER:
MATTHEWS/T/E
OFFICER I.D.: 582
PERMIT NUMBER: 220253
EXPIRATION DATE: 09/11/14
MISCELLANEOUS DATA:
RFI TEST

--- BREATH ANALYSIS ---

BLANK TEST	.000	00:00
INTERNAL STANDARD	VERIFIED	00:00
RADIO INTERFERENCE		

Operator Signature

Jim Matthews

2200-02

STATE OF MISSISSIPPI
DEPARTMENT OF HEALTH



PERMIT
TYPE II



TIMOTHY E MATTHEWS

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field repairs, and operate the following breath analyzer(s):

DATAMASTER

for the determination of the alcoholic content of blood from a sample of expired (alveolar) air. Issued under the provisions of sections 577.020 through 577.041, RSMo. 1986.

Date 09/11/2012

Number 220253

Expires 09/11/2014

Director of State Public Health Laboratory

Director, Department of Health



Simulator Calibration Report

This calibration report is to certify the alcohol reference simulator listed below has been examined and tested using standards traceable to the National Institute of Standards and Technology (NIST) in accordance with the standards set by the Missouri Department of Health and Senior Services Rules and Regulations:
19 CSR 25-30.051 (4).

Checked: 08/26/2014 Expires: 08/26/2015
Digital Therm. SN:094948
MSC Tech:RWV Temp:34.01
Agency: Hollister Police Dept
SD 2732



Technician Printed Name: ROBERT WELSH

Technician Signature: *Robert Welsh*

Date: 08/26/2014

Contact: Missouri Safety Center
Breath-Alcohol Instrument Training Program
660-543-4834