



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
 STATE PUBLIC HEALTH LABORATORY
 BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

RECEIVED
 By Carol Day at 12:21 pm, Oct 28, 2014
 REPORT #6

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).
 Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.
 Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 201229	NAME OF AGENCY St Robert Police Department	DATE OF INSPECTION 10/16/2014
LOCATION OF INSTRUMENT (STREET AND CITY) 194 Eastlawn Ave Ste A St Robert, MO 65584		TIME OF INSPECTION 4:26 pm

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) 10-16-2014/1626
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER _____ 49 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

INDICATOR LIGHTS

SIMULATOR SOLUTION SUPPLIER Guth Laboratories LOT # 13100 EXP. DATE 04/23/2015

SIMULATOR TEMP (34°C ± 0.2°C) 34.0 °C SIMULATOR SN SD2284 EXP. DATE 07/15/2015

CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

- 0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
- 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
- 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 <input checked="" type="checkbox"/> 0.097%	TEST 2 <input checked="" type="checkbox"/> 0.098%	TEST 3 <input checked="" type="checkbox"/> 0.098%
---	---	---

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:
 (DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	0	(0-.04)	0	(.05-.09)	1	(.10-.14)	4	(.15-.19)	2	OVER .19	0
----------	---	---------	---	-----------	---	-----------	---	-----------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS
 (USE OTHER SIDE IF NECESSARY).

Instrument is operating within requirements of Missouri DHSS

INSPECTING OFFICER

SIGNATURE 	PRINT FULL NAME Jennifer R Janko
TYPE II PERMIT NUMBER/EXPIRATION DATE 240037 02/24/2016	TELEPHONE NUMBER (573) 336-4700

RETURN COMPLETED REPORT TO THE:
 Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
 2875 James Blvd.
 Poplar Bluff, MO 63901



GUTH LABORATORIES, INC.

590 NORTH 67th STREET • HARRISBURG, PA 17111-4611 • TELEPHONE: 717-564-6470

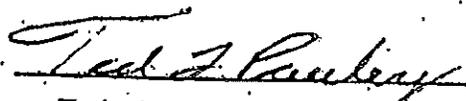
CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 13100 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on April 29, 2013, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1214% (w/vol) ethyl alcohol. The expiration date for this lot number is April 23, 2015 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.


Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Certified Reference Standard lot number FN122211-02 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.



STATE OF MISSOURI
DEPARTMENT OF HEALTH AND SENIOR SERVICES
BREATH ALCOHOL PROGRAM



**PERMIT
TYPE II**

JENNIFER R JANKO

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

DATAMASTER

for the determination of the alcoholic content of blood from a sample of expired air; Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE: 2/24/2014

NUMBER 240037

EXPIRES 2/24/2016

MO-530-0771 (6/10)

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

LAB-4 (R6-10)



STATE OF MISSOURI
DEPARTMENT OF HEALTH AND SENIOR SERVICES
BREATH ALCOHOL PROGRAM

INSTRUMENT OPERATOR CARD

The named cardholder is authorized to operate an evidential breath alcohol instrument for the determination of the alcoholic content in breath form of expired air in Missouri.



Operator: JANKO, JENNIFER
Permit No. 240037
Date Issued 2/24/2014 Date Expires 2/24/2016

BAC DataMaster
Evidence Ticket

STATE OF MISSOURI

MOO DRAINAGE DISTRICT SERIAL NUMBER 201229
10/16/12

ARREST TIME: 00:00

SUBJECT NAME:

RF CHECK

DOB: 01/01/01

SEX: F

STATE/M.L.: MO/NONE

ARRESTING OFFICER:

MR

OFFICER I.D.: 000

TESTING OFFICER:

JANIKO/JENNIFER/R

OFFICER I.D.: 100

PERMIT NUMBER: 240007

EXPIRATION DATE: 02/24/16

RESOLUTIONS DATA:

--- BREATH ANALYSIS ---

BUNK TEST

.000

16:37

INTERNAL STANDARD

VERIFIED

16:37

AUDIO INTERFERENCE

Operator Signature

