



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES  
 STATE PUBLIC HEALTH LABORATORY  
 BREATH ALCOHOL PROGRAM  
**DATAMASTER MAINTENANCE REPORT**

**RECEIVED**  
 By Carol Day at 12:04 pm, Dec 23, 2014  
 REPORT #6

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).  
 Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.  
 Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 201248	NAME OF AGENCY Cameron P.D.	DATE OF INSPECTION 12-17-2014
LOCATION OF INSTRUMENT (STREET AND CITY) 101 North Chestnut, Cameron, Missouri, 64429		TIME OF INSPECTION 10:55

**CHECKLIST:** Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) 12-17-2014 10:55
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 49°C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

<input checked="" type="checkbox"/> INDICATOR LIGHTS
<input checked="" type="checkbox"/> SIMULATOR SOLUTION SUPPLIER Guth Laboratories Inc. LOT # 13290 EXP. DATE 10-29-2015
<input checked="" type="checkbox"/> SIMULATOR TEMP (34°C ± 0.2°C) 34 °C SIMULATOR SN SD2263 EXP. DATE 08/01/2015

CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE  
 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE  
 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 .095	TEST 2 .096	TEST 3 .096
-------------	-------------	-------------

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:  
 (DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	0	(0-.04)	0	(.05-.09)	0	(.10-.14)	0	(.15-.19)	0	OVER .19	0
----------	---	---------	---	-----------	---	-----------	---	-----------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

Tested and certified with in DOH requirements

INSPECTING OFFICER	
SIGNATURE 	PRINT FULL NAME Corporal Shane T. Stevenson
TYPE II PERMIT NUMBER/EXPIRATION DATE 230138/07-17-15	TELEPHONE NUMBER (816) 632-6521

RETURN COMPLETED REPORT TO THE: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office  
 2875 James Blvd.  
 Poplar Bluff, MO 63901



**GUTH LABORATORIES, INC.**

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-634-5470.

## CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 13290 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on October 31, 2013, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1202% (w/vol) ethyl alcohol. The expiration date for this lot number is October 29, 2015 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President  
GUTH LABORATORIES, INC.

**NIST Traceability:**

Testing was conducted using Cerilliant Reference Standard lot number FN122211-02 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.



STATE OF MISSOURI  
DEPARTMENT OF HEALTH AND SENIOR SERVICES  
BREATH ALCOHOL PROGRAM

2

PERMIT  
TYPE II

SHANE T STEVENSON

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

DATAMASTER

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 7/17/2013

NUMBER 230138

EXPIRES 7/17/2015

MO 530-7771 (8-10)

*[Signature]*

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

*[Signature]*

acting director

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES


LAB-4 (RS-10)



STATE OF MISSOURI  
DEPARTMENT OF HEALTH AND SENIOR SERVICES  
BREATH ALCOHOL PROGRAM

INSTRUMENT OPERATOR CARD

The named cardholder is authorized to operate an evidential breath alcohol instrument for the determination of the alcoholic content in breath form of expired air in Missouri.

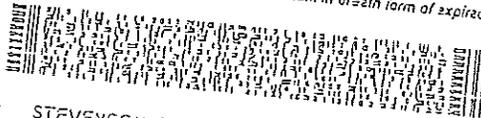

Operator STEVENSON, SHANE  
Permit No 230138  
Date Issued 7/17/2013 Date Expires 7/17/2015

Face This Side Down - This Edge In First

**BAC DataMaster**  
Evidence Ticket

Face This Side Down - This Edge In First

**BAC DataMaster**  
Evidence Ticket

Signature \_\_\_\_\_

Operator Signature \_\_\_\_\_

Face This Side Down - This Edge In First

**BAC DataMaster**  
Evidence Ticket

Operator Signature \_\_\_\_\_

*[Handwritten Signature]*