



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
 STATE PUBLIC HEALTH LABORATORY
 BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

RECEIVED REPORT #6
 By Carol Day at 9:45 am, Aug 28, 2014

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).
 Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.
 Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 201294	NAME OF AGENCY MALDEN POLICE DEPARTMENT	DATE OF INSPECTION 08-20-2014
LOCATION OF INSTRUMENT (STREET AND CITY) 301 E. BOURKE ST. MALDEN MD 63552		TIME OF INSPECTION 11:12

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) 08/20/14 11:12
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 49 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

INDICATOR LIGHTS

SIMULATOR SOLUTION SUPPLIER **LUTIX LAB** LOT # **13290** EXP. DATE **10-29-15**

SIMULATOR TEMP (34°C ± 0.2°C) **34°** °C SIMULATOR SN **S02666B** EXP. DATE **2-14-15**

CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 • .099	TEST 2 • .101	TEST 3 • .101
----------------------	----------------------	----------------------

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:
 (DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS <input checked="" type="checkbox"/> (0-.04) 1	(.05-.09) <input checked="" type="checkbox"/>	(.10-.14) 2	(.15-.19) <input checked="" type="checkbox"/>	OVER .19 1
---	---	--------------------	---	-------------------

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

MEETS DOH STANDARDS

INSPECTING OFFICER

SIGNATURE Jeff Brown	PRINT FULL NAME JEFF BROWN
TYPE II PERMIT NUMBER/EXPIRATION DATE 240263 6-2-16	TELEPHONE NUMBER 660 385 2195

RETURN COMPLETED REPORT TO THE: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
 2875 James Blvd.
 Poplar Bluff, MO 63901



GUTH LABORATORIES, INC.

889 NORTH 67th STREET • HARRISBURG, PA. 17111-4611 • TELEPHONE: 717-664-5470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 13290 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on October 31, 2013, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1202% (w/vol) ethyl alcohol. The expiration date for this lot number is October 29, 2015 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Certified Reference Standard lot number FN122211-02 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.



STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

2

**PERMIT
 TYPE II**

JEFF BROWN

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

DATAMASTER

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 6/2/2014

NUMBER 240263

EXPIRES 6/2/2016

MO.590-9771 (6-10)

W. S. ...

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

Paul Vackerly

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

LAB-4 (06-10)

FACE THIS SIDE DOWN - THIS EDGE IN FIRST

BAC DataMaster
Evidence Ticket

STATE OF MISSOURI
MALDEN POLICE DEPARTMENT

BAC DATA MASTER SERIAL NUMBER 201294
08/26/14

ARREST TIME: 11:09

SUBJECT NAME:

WLOW/DE/S

DUB: 12/04/08 SEX: M

STATE/I.L.: MO/M495800782

ARRESTING OFFICER:

BROWN/JEFF

OFFICER I.D. # 765 7 11 11

TESTING OFFICER:

S/A/R

OFFICER I.D.: S/A/R

PERMIT NUMBER: 240563

EXPIRATION DATE: 06/02/16

MISCELLANEOUS DATA:


RFI CHECK

N/A

--- BREATH ANALYSIS ---

BLANK TEST	1.000	11:21
INTERNAL STANDARD	VERIFIED	11:22
RADIO INTERFERENCE		

OPERATOR SIGNATURE

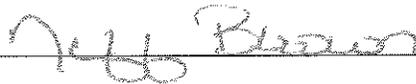

Card Stock No.
60021

REORDER ALL SUPPLIES FROM N.P.A.S.
P.O. BOX 1435, MANSFIELD, OH 44901