



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
 STATE PUBLIC HEALTH LABORATORY
 BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

RECEIVED
 By Carol Day at 1:13 pm, Oct 16, 2014

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).
 Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.
 Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 202049	NAME OF AGENCY Waynesville Police Department	DATE OF INSPECTION 10/15/2014
LOCATION OF INSTRUMENT (STREET AND CITY) 201 North Street Waynesville		TIME OF INSPECTION 6:11 pm

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) <u>10/15/2014 17:41</u>
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER <u>49</u> °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

INDICATOR LIGHTS

SIMULATOR SOLUTION SUPPLIER RepCo Marketing Inc. LOT # 13002 EXP. DATE 06/19/2015

SIMULATOR TEMP (34°C ± 0.2°C) 34 °C SIMULATOR SN SD2283 EXP. DATE 02/10/2015

CALIBRATION CHECK – (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

- 0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
- 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
- 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 .097	TEST 2 .098	TEST 3 .098
--------------	--------------	--------------

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:
 (DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	0	(.0-.04)	6	(.05-.09)	1	(.10-.14)	0	(.15-.19)	0	OVER .19	0
----------	---	----------	---	-----------	---	-----------	---	-----------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

INSPECTING OFFICER

SIGNATURE 	PRINT FULL NAME Corporal John Meir
TYPE (PERMIT NUMBER/EXPIRATION DATE) 240343 07/30/2016	TELEPHONE NUMBER (573) 774-2414

RETURN COMPLETED REPORT TO THE: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
 2875 James Blvd.
 Poplar Bluff, MO 63901

CERTIFICATE OF ANALYSIS

MANUFACTURER AND SUPPLIER: RepCo Marketing, Inc.

LOT NUMBER: 13002

EXPIRATION DATE: June 19, 2015 at 11:59 p.m.

RepCo Marketing, Inc. certifies the following:


RepCo Marketing, Inc. manufactured, tested and supplied Lot Number 13002 of Alcohol Certified Solution for simulators. Random samples of said lot number were analyzed by an independent laboratory utilizing a gas chromatograph and found to contain .1217 gms/dl +/- .003 gms/dl wt/vol ethanol (95% Confidence).

The alcohol and distilled water used in the solution were found to be free of any interfering substance.

This solution will produce a vapor alcohol value of .100 +/- 3% gms/210L Breath when heated to 34 Degrees Celsius +/- 0.2 Degrees Celsius in a simulator (95% Confidence).

The date of manufacture for this lot number is June 20, 2013
The expiration date for this lot number is June 19, 2015 at
11:59 p.m.

This document is a true representation of the original Certificate of Analysis.

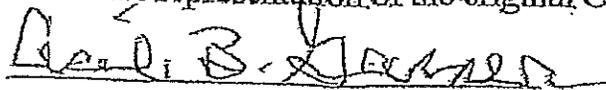

Cecil B. Garner, President
RepCo Marketing, Inc.

**BAC DataMaster
Evidence Ticket**

STATE OF MISSOURI
WAYNESVILLE POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 202049
10/15/14

ARREST TIME: 17:15
SUBJECT NAME:
RFI/TEST
DOB: 10/15/14 SEX: M
STATE/D.L.: MO/12345
ARRESTING OFFICER:
MEIR/J
OFFICER I.D.: 705
TESTING OFFICER:
MEIR/J
OFFICER I.D.: 705
PERMIT NUMBER: 240313
EXPIRATION DATE: 07/30/16
MISCELLANEOUS DATA:

--- BREATH ANALYSIS ---

BLANK TEST .000 17:45
INTERNAL STANDARD VERIFIED 17:45
RADIO INTERFERENCE

Operator Signature



2200-02

**BAC DataMaster
Evidence Ticket**

STATE OF MISSOURI
WAYNESVILLE POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 202049
10/15/14

TESTING OFFICER:
MEIR/JOHN
OFFICER I.D.: 705
PERMIT NUMBER: 240313
EXPIRATION DATE: 07/30/16
MISCELLANEOUS DATA:

--- SUPERVISOR MODE ---

BLANK TEST .000 18:08
INTERNAL STANDARD VERIFIED 18:08
EXTERNAL STANDARD .097 18:08
BLANK TEST .000 18:09
EXTERNAL STANDARD .098 18:09
BLANK TEST .000 18:10
EXTERNAL STANDARD .098 18:11
BLANK TEST .000 18:11

N = 3
SIM. = .1
AVG. = .0976

Operator Signature



2208

Face This Side Down - This Edge In First

BAC DataMaster
Evidence Ticket

STATE OF MISSOURI
WAYNESVILLE POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 202049
10/15/14
17:41

--- DIAGNOSTIC CHECK ---

COMPUTER: OKAY

PROGRAM (04-07-2009): OKAY

HEATERS
SAMPLE CHAMBER: 49c

FLOW DETECTOR: OKAY

PUMP
HIGH SPEED: OKAY

DETECTOR: OKAY

FILTERS: OKAY


QUARTZ STANDARD: OKAY

CALIBRATION: OKAY

PRINTER TEST

!"#\$%&'()*+,-./0123456789:;<=>?@ABCDEFGH
IJKLMNOPQRSTUVWXYZ[\]^_`abode fghijk lmn
opqrstuvwxyz{|}~

Operator Signature

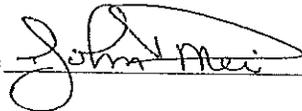



STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

2

PERMIT
TYPE II
JOHN P MEIR

Is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

DATAMASTER, INTOX EC/IR II

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 7/30/2014

NUMBER 240313

EXPIRES 7/30/2016

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

MO 680-0771 (0-10)

LAB-4 (R3-10)

STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

INSTRUMENT OPERATOR CARD

The named cardholder is authorized to operate an evidential breath alcohol instrument for the determination of the alcoholic content in breath form of expired air in Missouri.

Operator MEIR, JOHN
 Permit No 240313
 Date Issued 7/30/2014 Date Expires 7/30/2016



Simulator Calibration Report

This is to certify that the Simulator listed below has been examined and tested using standards traceable to the National Institute of Standards and Technology in accordance to the standards set by Missouri Department of Health and Senior Services Rules and Regulations 19CSR 25-30

Checked: 2/10/2014 Expires: 02/10/2015
Digital Therm. SN:094948 Temp:33.96
MSC Tech:DRL
Agency: Waynesville PD



SD2283



Tested By:

DAN LUCAS

Date:

2/10/2014