



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
 STATE PUBLIC HEALTH LABORATORY
 BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

REPORT #6

RECEIVED
 By Carol Day at 11:36 am, Dec 18, 2014

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 30 days) or whenever the instrument is serviced or repaired and whenever it is placed into service. Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 204112	NAME OF AGENCY MSHP	DATE OF INSPECTION 12/04/2014
LOCATION OF INSTRUMENT (STREET AND CITY) RANDOLPH CO SO 372 HWY JJ HUNTSVILLE		TIME OF INSPECTION 2:45 pm

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) 12/04/2014 1444
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 50 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

<input checked="" type="checkbox"/> INDICATOR LIGHTS
<input checked="" type="checkbox"/> SIMULATOR SOLUTION SUPPLIER Guth LOT # 13290 EXP. DATE 10/29/2015
<input checked="" type="checkbox"/> SIMULATOR TEMP (34°C ± 0.2°C) 34.0 °C SIMULATOR SN MP2518 EXP. DATE 09/19/2015

CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

<input checked="" type="checkbox"/> 0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
<input type="checkbox"/> 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
<input type="checkbox"/> 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 .097%	TEST 2 .096%	TEST 3 .096%
--------------	--------------	--------------

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT: (DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	0	(0-.04)	0	(.05-.09)	2	(.10-.14)	1	(.15-.19)	0	OVER .19	0
----------	---	---------	---	-----------	---	-----------	---	-----------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

INSPECTING OFFICER	
SIGNATURE 	PRINT FULL NAME Bruce E. Thoman
TYPE II PERMIT NUMBER/EXPIRATION DATE 240185 04/22/2016	TELEPHONE NUMBER (660) 385-2132

RETURN COMPLETED REPORT TO THE: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
 2875 James Blvd.
 Poplar Bluff, MO 63901



GUTH LABORATORIES, INC.

660 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 13290 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on October 31, 2013, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1202% (w/vol) ethyl alcohol. The expiration date for this lot number is October 29, 2015 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FNI22211-02 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

BAC DataMaster
Evidence Ticket

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204112
12/04/14

ARREST TIME: 10:00
SUBJECT NAME:
RFI/TEST
DOB: 11/11/11 SEX: M
STATE/D.L.: MO/111111
ARRESTING OFFICER:
*
OFFICER I.D.: *
TESTING OFFICER:
THOMAS/BRUCE/E
OFFICER I.D.: 1256
PERMIT NUMBER: 240185
EXPIRATION DATE: 04/22/16
MISCELLANEOUS DATA:

--- BREATH ANALYSIS ---

BLANK TEST	.000	15:00
INTERNAL STANDARD	VERIFIED	15:00
RADIO INTERFERENCE		



Operator Signature _____

BAC DataMaster
Evidence Ticket

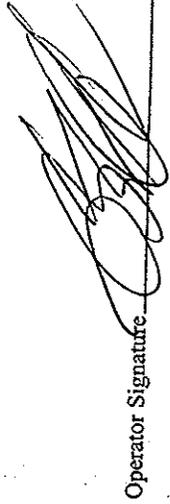
MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204112
12/04/14

TESTING OFFICER:
THOMAS/BRUCE/E
OFFICER I.D.: 1256
PERMIT NUMBER: 240185
EXPIRATION DATE: 04/22/16
MISCELLANEOUS DATA:

--- SUPERVISOR MODE ---

BLANK TEST	.000	14:54
INTERNAL STANDARD	VERIFIED	14:54
EXTERNAL STANDARD	.097	14:55
BLANK TEST	.000	14:55
EXTERNAL STANDARD	.096	14:56
BLANK TEST	.000	14:56
EXTERNAL STANDARD	.096	14:57
BLANK TEST	.000	14:57

N = 3
SIM. = .1
AVG. = .0963

Operator Signature 

BAC DataMaster
Evidence Ticket

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204112
12/04/14
14:44

--- DIAGNOSTIC CHECK ---

COMPUTER: OKAY

PROGRAM (04-07-2009): OKAY

HEATERS

SAMPLE CHAMBER: 50c

FLOW DETECTOR: OKAY

PUMP

HIGH SPEED: OKAY

DETECTOR: OKAY

FILTERS: OKAY

QUARTZ STANDARD: OKAY

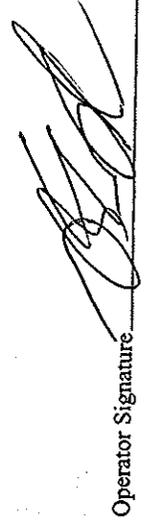
CALIBRATION: OKAY

PRINTER TEST

! " # \$ % & ' () * + , - . / 0 1 2 3 4 5 6 7 8 9 ; : < = > ? @ A B C D E F G

H I J K L M N O P Q R S T U V W X Y Z [\] ^ _ ` a b c d e f g h i j k l m n o

p q r s t u v w x y z { | } ~

Operator Signature 



STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

2

PERMIT
TYPE II
BRUCE E THOMAN

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

DATAMASTER, INTOX DMT

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 4/22/2014

NUMBER 240185

EXPIRES 4/22/2016

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

MO 580-0771 (6-10)

LAB-4 (R6-10)

STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

INSTRUMENT OPERATOR CARD

The named cardholder is authorized to operate an evidential breath alcohol instrument for the determination of the alcoholic content in breath form of expired air in Missouri.

Operator THOMAN, BRUCE
 Permit No 240185
 Date Issued 4/22/2014 Date Expires 4/22/2016