



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
 STATE PUBLIC HEALTH LABORATORY
 BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

received 9/15/14-cd **REPORT #6**

REVIEWED
 By Carol Day at 11:27 am, Sep 25, 2014

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 30 days). Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service. Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 204129	NAME OF AGENCY Missouri State Highway Patrol	DATE OF INSPECTION 09/07/2014
LOCATION OF INSTRUMENT (STREET AND CITY) 891 Technology Dr., Weldon Spring, MO - Troop C "Road Room"		TIME OF INSPECTION 11:05 am

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) <u>09/07/2014 11:05 am</u>
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER <u>49 °C</u>	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER
<input checked="" type="checkbox"/> INDICATOR LIGHTS	
<input checked="" type="checkbox"/> SIMULATOR SOLUTION SUPPLIER <u>Guth Laboratories, Inc</u>	LOT # <u>13280</u> EXP. DATE <u>10/16/2015</u>
<input checked="" type="checkbox"/> SIMULATOR TEMP (34°C ± 0.2°C) <u>34.0</u> °C	SIMULATOR SN <u>G11040</u> EXP. DATE <u>05/22/2015</u>

CALIBRATION CHECK – (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 <input checked="" type="checkbox"/> .099%	TEST 2 <input checked="" type="checkbox"/> .099%	TEST 3 <input checked="" type="checkbox"/> .100%
--	--	--

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT: (DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	0	(0-.04)	0	(.05-.09)	0	(.10-.14)	0	(.15-.19)	0	OVER .19	0
----------	---	---------	---	-----------	---	-----------	---	-----------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

This unit meets all Department of Health rules and regulations.
 Bottle #643 Inventory #127255

INSPECTING OFFICER		PRINT FULL NAME J.T. Hedrick
SIGNATURE 	TELEPHONE NUMBER (636) 300-2800	
TYPE II PERMIT NUMBER/EXPIRATION DATE 240055	03/07/2016	

RETURN COMPLETED REPORT TO THE: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
 2875 James Blvd.
 Poplar Bluff, MO 63901

Face This Side Down - This Edge In First

**BAC DataMaster
Evidence Ticket**

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204129
09/07/14
11:05

--- DIAGNOSTIC CHECK ---

COMPUTER: OKAY
PROGRAM (04-07-2009): OKAY
HEATERS
SAMPLE CHAMBER: 49c
FLOW DETECTOR: OKAY
PUMP
HIGH SPEED: OKAY
DETECTOR: OKAY
FILTERS: OKAY
QUARTZ STANDARD: OKAY
CALIBRATION: OKAY

PRINTER TEST

!"#\$%&'()*+,-./0123456789:;<=>?@ABCDEFGH
IJKLMNOPQRSTUVWXYZ[\]^_`abcdefghijklmnopq
rstuvwxyz{|}~

Operator Signature



Face This Side Down - This Edge In First

**BAC DataMaster
Evidence Ticket**

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204129
09/07/14


TESTING OFFICER:
HEDRICK/J/T
OFFICER I.D.: 455
PERMIT NUMBER: 240055
EXPIRATION DATE: 03/07/16
MISCELLANEOUS DATA:

--- SUPERVISOR MODE ---

BLANK TEST	.000	11:11
INTERNAL STANDARD	VERIFIED	11:11
EXTERNAL STANDARD	.099	11:11
BLANK TEST	.000	11:12
EXTERNAL STANDARD	.099	11:12
BLANK TEST	.000	11:13
EXTERNAL STANDARD	.100	11:13
BLANK TEST	.000	11:14

N = 3
SIM. = .1
AVG. = .0993

Operator Signature

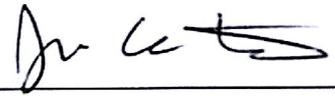

Face This Side Down - This Edge In First

BAC DataMaster
Evidence Ticket


MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204129
09/07/14

ARREST TIME: 09:00
SUBJECT NAME:
TEST
DOB: 11/11/11 SEX: M
STATE/D.L.: MO/1234567890
ARRESTING OFFICER:
TEST
OFFICER I.D.: 1
TESTING OFFICER:
HEDRICK/J/T
OFFICER I.D.: 455
PERMIT NUMBER: 240055
EXPIRATION DATE: 03/07/16
MISCELLANEOUS DATA:
RFI TEST

--- BREATH ANALYSIS ---

BLANK TEST	.000	11:17
INTERNAL STANDARD	VERIFIED	11:17
RADIO INTERFERENCE		

Operator Signature _____

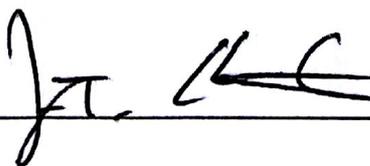



STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

2

**PERMIT
 TYPE II**

JAMES T HEDRICK

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

DATAMASTER, INTOX DMT

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 3/7/2014

NUMBER 240055

EXPIRES 3/7/2016

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

LAB-4 (R6-10)

MO 580-0771 (6-10)



GUTH LABORATORIES, INC.

800 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 13280 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on October 18, 2013, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1217% (w/vol) ethyl alcohol. The expiration date for this lot number is October 16, 2015 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN122211-02 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.