



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
STATE PUBLIC HEALTH LABORATORY
BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

REPORT #6

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days). Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service. Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 204190 Inv# 127310	NAME OF AGENCY MSHP	DATE OF INSPECTION 04/01/2014
LOCATION OF INSTRUMENT (STREET AND CITY) Polk County Jail, Bolivar		TIME OF INSPECTION 10:16 am

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) 04/01/2014 10:16
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER _____ 49 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

INDICATOR LIGHTS

SIMULATOR SOLUTION SUPPLIER GUTH LABORATORIES INC. LOT # 13280 EXP. DATE 10/16/2015

SIMULATOR TEMP (34°C ± 0.2°C) 34.15 °C SIMULATOR SN G6804 EXP. DATE 07/15/2014

CALIBRATION CHECK – (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

- 0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
- 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
- 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 .099	TEST 2 .099	TEST 3 .100
--------------	--------------	--------------

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT: (DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	0	(0-.04)	1	(.05-.09)	3	(.10-.14)	3	(.15-.19)	1	OVER .19	0
----------	---	---------	---	-----------	---	-----------	---	-----------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

INSPECTING OFFICER

SIGNATURE 	PRINT FULL NAME Joshua L. White
---------------	------------------------------------

TYPE II PERMIT NUMBER/EXPIRATION DATE 240080 03/07/2016	TELEPHONE NUMBER (417) 895-6868
--	------------------------------------

RETURN COMPLETED REPORT TO THE: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
2875 James Blvd.
Poplar Bluff, MO 63901



GUTH LABORATORIES, INC.

880 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-664-6470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 13280 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on October 18, 2013, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1217% (w/vol) ethyl alcohol. The expiration date for this lot number is October 16, 2015 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FNI22211-02 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

Face This Side Down - This Edge In First

**BAC DataMaster
Evidence Ticket**

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204190
04/01/14

TESTING OFFICER:
WHITE/JOSHUA/L
OFFICER I.D.: 1202
PERMIT NUMBER: 240030
EXPIRATION DATE: 03/07/16
MISCELLANEOUS DATA:
MAINTENANCE

--- SUPERVISOR MODE ---

BLANK TEST	.000	10:20
INTERNAL STANDARD	VERIFIED	10:20
EXTERNAL STANDARD	.099	10:20
BLANK TEST	.000	10:21
EXTERNAL STANDARD	.099	10:22
BLANK TEST	.000	10:22
EXTERNAL STANDARD	.100	10:23
BLANK TEST	.000	10:23

N = 3
SIM. = .1
AVG. = .0993

Operator Signature



Face This Side Down - This Edge In First

**BAC DataMaster
Evidence Ticket**

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204190
04/01/14
10:16


--- DIAGNOSTIC CHECK ---

COMPUTER:	OKAY
PROGRAM (04-07-2009):	OKAY
HEATERS	
SAMPLE CHAMBER:	49c
FLOW DETECTOR:	OKAY
PUMP	
HIGH SPEED:	OKAY
DETECTOR:	OKAY
FILTERS:	OKAY
QUARTZ STANDARD:	OKAY
CALIBRATION:	OKAY

PRINTER TEST

!"#\$%&'()*+,-./0123456789:;<=>?@ABCDEFGHIJ
KLMNOPQRSTUVWXYZ[\]^_`abcdefghijklmnopq
rstuvwxyz{|}~

Operator Signature

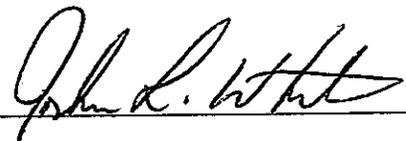

Face This Side Down - This Edge in First

**BAC DataMaster
Evidence Ticket**


MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204190
04/01/14

ARREST TIME: 01:23
SUBJECT NAME:
X/X/X
DOB: 01/23/45 SEX: M
STATE/D.L.: MO/XXXXXXXXXX
ARRESTING OFFICER:
SAME
OFFICER I.D.: SAME
TESTING OFFICER:
WHITE/JOSHUA/L
OFFICER I.D.: 1202
PERMIT NUMBER: 240080
EXPIRATION DATE: 03/07/16
MISCELLANEOUS DATA:
MAINTENANCE

--- BREATH ANALYSIS ---

BLANK TEST	.000	10:26
INTERNAL STANDARD	VERIFIED	10:27
RADIO INTERFERENCE		

Operator Signature

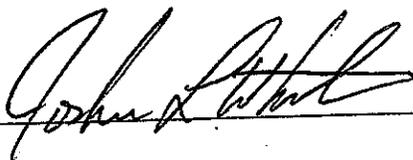



STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

2

**PERMIT
 TYPE II**

JOSHUA L WHITE

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

DATAMASTER, INTOX DMT, ALCO-SENSOR IV W/PRINTER

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 3/7/2014

NUMBER 240080

EXPIRES 3/7/2016

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

MO 580-0771 (6-10)

LAB-4 (R6-10)



STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

INSTRUMENT OPERATOR CARD

The named cardholder is authorized to operate an evidential breath alcohol instrument for the determination of the alcoholic content in breath form of expired air in Missouri.



Operator **WHITE, JOSHUA**
 Permit No **240080**
 Date Issued **3/7/2014** Date Expires **3/7/2016**