



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
 STATE PUBLIC HEALTH LABORATORY
 BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

DOH

RECEIVED
 By Carol Day at 2:37 pm, Oct 17, 2014 RT #6

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).
 Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.
 Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 204191	NAME OF AGENCY Missouri State Highway Patrol	DATE OF INSPECTION 10/02/2014
LOCATION OF INSTRUMENT (STREET AND CITY) Newton County Sheriffs Department, Neosho		TIME OF INSPECTION 2:20 am

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) 10/02/2014 @ 0224 hours
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 50°C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

INDICATOR LIGHTS

SIMULATOR SOLUTION SUPPLIER Guth Laboratories LOT # 13290 EXP. DATE 10/29/2015

SIMULATOR TEMP (34°C ± 0.2°C) 34.0 °C SIMULATOR SN SD3143 EXP. DATE 01/16/2015

CALIBRATION CHECK – (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

- 0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
- 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
- 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 ← .101	TEST 2 ← .102	TEST 3 ← .102
---------------	---------------	---------------

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:
 (DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	3	(0-.04)	0	(.05-.09)	1	(.10-.14)	3	(.15-.19)	7	OVER .19	2
----------	---	---------	---	-----------	---	-----------	---	-----------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

INSPECTING OFFICER	
SIGNATURE 	PRINT FULL NAME Jeffrey T Day
TYPE II PERMIT NUMBER/EXPIRATION DATE 230019 / 02/08/2015	TELEPHONE NUMBER (417) 895-6868

RETURN COMPLETED REPORT TO THE:
 Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
 2875 James Blvd.
 Poplar Bluff, MO 63901



GUTH LABORATORIES, INC.

580 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-684-5470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 13290 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on October 31, 2013, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1202% (w/vol) ethyl alcohol. The expiration date for this lot number is October 29, 2015 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.


The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN122211-02 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

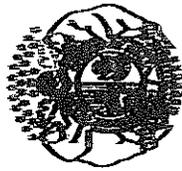

STATE OF MISSOURI
DEPARTMENT OF HEALTH AND SENIOR SERVICES
BREATH ALCOHOL PROGRAM

2

PERMIT
TYPE II

JEFFREY T DAY

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

DATAMASTER

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 308.111 through 308.113 RSMo.

DATE 02/08/2013

NUMBER 230019

EXPIRES 02/08/2015

W. W. S.

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

Richard W. S. Vandenberg

Acting Director

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

MS 610-0771 (8-10)

LAB-6 (HR-10)

Face This Side Down - This Edge In First

**BAC DataMaster
Evidence Ticket**

MISSOURI STATE HIGHWAY FINROU
BAC DATAMASTER SERIAL NUMBER 044191
10/09/14

TECHNIC OFFICER:
DAY/MT
OFFICER I.D.# 360
PERMIT NUMBER: 224415
EXPIRATION DATE: 09/08/15
MISCELLANEOUS DATA:

--- SUPERVISOR CODE ---

ALUMINUM TEST	1000	00100
INTERNAL STANDARD	954110	00100
EXTERNAL STANDARD	1191	00100
ALUMINUM TEST	1000	00100
INTERNAL STANDARD	1100	00100
ALUMINUM TEST	000	00100
EXTERNAL STANDARD	1100	00100
ALUMINUM TEST	000	00100

W = 0
STL = 0
RYS = 1100

Operator Signature

