



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES  
STATE PUBLIC HEALTH LABORATORY  
BREATH ALCOHOL PROGRAM

INTOX DMT MAINTENANCE REPORT

RECEIVED  
By Carol Day at 11:50 am, May 07, 2014

REPORT #1

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).  
Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.  
Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

INTOX DMT SN 500099	NAME OF AGENCY Missouri State Highway Patrol	DATE OF INSPECTION 04/29/2014
LOCATION OF INSTRUMENT (STREET AND CITY) 1510 E. Elm Street, Jefferson City		TIME OF INSPECTION 14:51:29

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or is operating within established limits. (Write in observed values where determined). Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC RECORD	
DATE AND TIME <u>04/29/2014 14:51:31</u>	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTER 1
<input checked="" type="checkbox"/> SAMPLE CHAMBER <u>48.9°C</u>	<input checked="" type="checkbox"/> FILTER 2
<input checked="" type="checkbox"/> BREATH TUBE <u>48.1°C</u>	<input checked="" type="checkbox"/> FILTER 3
<input checked="" type="checkbox"/> PUMP	<input checked="" type="checkbox"/> INTERNAL STANDARD

BREATH ANALYZER ACCURACY STANDARDS

<input type="checkbox"/> SIMULATOR STANDARD	<input checked="" type="checkbox"/> COMPRESSED ETHANOL-GAS MIXTURE
<input checked="" type="checkbox"/> STANDARD SUPPLIER <u>ILMO</u>	LOT # <u>21913080A4</u> EXP. DATE <u>09/01/2015</u>
<input type="checkbox"/> SIMULATOR TEMP (34°C ± 0.2°C) _____	SIMULATOR SN _____ SIMULATOR EXP DATE _____

CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)  
Run three tests using a standard. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard being used.

<input type="checkbox"/> 0.10% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
<input checked="" type="checkbox"/> 0.08% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
<input type="checkbox"/> 0.04% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1: 0.076	TEST 2: 0.077	TEST 3: 0.077
---------------	---------------	---------------

PERFORM R.F.I. TEST

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:

REFUSALS: 1	0-.04: 1	.05-.09: 0	.10-.14: 4	.15-.19: 0	OVER .19: 2
-------------	----------	------------	------------	------------	-------------

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY)

---



---



---



---

INSPECTING OFFICER

SIGNATURE 	PRINT FULL NAME MARK A WINDER	
TYPE II PERMIT NUMBER 240082	EXPIRATION DATE 03/07/2016	TELEPHONE NUMBER

RETURN COMPLETED REPORT TO THE Breath Alcohol Program, MO Department of Health and Senior Services  
Southeast District Office  
2875 James Blvd, Poplar Bluff, MO 63901



7 Eastgate Dr. • P.O. Box 790 • Jacksonville, IL 62651-0790  
217-245-2183 • Fax: 217-243-7634 • www.ilmoproducts.com

### Certificate of Analysis

Certificate ID: 5180  
Part #: BAC105L080T  
Cylinder Size: 105L  
Lot Number: 21913080A4  
Expiration: 9/1/2015

0.080 BAC (For use with breath alcohol testing instruments)


Contents: 105 Liters @ 1000 psig 70°F (21°C)

Component:	Concentration:	Accuracy:	Method:
Ethanol	208.4 ppm	+/- 0.002 or 2%	NDIR
Nitrogen	Balance	BAC whichever is greater	

\*NIST Standard Reference Material  
Cylinder No. CC14290 / Job No. 09160202  
Certified 212.8 µmol/mol Ethanol in Nitrogen  
for ILMO Products Co., Jacksonville, IL

*James Mattes*  
Specialty Gas Lab Tech

08/22/13  
Date

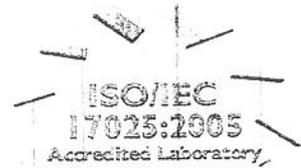

Distributed by: CMI Inc.  
316 East Ninth Street  
Owensboro, KY 42303  
Phone 866-835-0690  
[www.alcoholtest.com](http://www.alcoholtest.com)