

32



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
 STATE PUBLIC HEALTH LABORATORY
CMI INTOXILYZER 8000 MAINTENANCE REPORT

RECEIVED
 By Carol Day at 11:52 am, Sep 10, 2014

REPORT #2

Complete this report in duplicate at the time of the regular monthly preventive maintenance check, and whenever instrument is repaired. Send one copy to Department of Health and Senior Services, and retain one copy in department file.

INSTRUMENT SERIAL NUMBER 80-005840	LOCATION OF INSTRUMENT KCMO POLICE DEPT.	DATE OF INSPECTION 09/02/2014	TIME OF INSPECTION 21:57
---------------------------------------	---	----------------------------------	-----------------------------

CALIBRATION CHECK RESULTS			CALIBRATION CHECK SUMMARY		
Test	g/210L	Time	STANDARD TYPE	STANDARD LOT #	STANDARD EXPIRATION DATE
Air Blank	0.000	21:58	DRY	05514080	04/01/2016
Cal Check	0.079	21:59	SIM TEMPERATURE	SIM SERIAL NUMBER	SIM CERTIFICATE EXPIRATION
Air Blank	0.000	21:59	N/A	N/A	N/A
Cal Check	0.079	22:00	STANDARD VALUE	STANDARD SUPPLIER	
Air Blank	0.000	22:00	0.080	CMI	
Cal Check	0.079	22:00	CALIBRATION CHECK RESULT 1		
Air Blank	0.000	22:00	0.079		
Cal Check	0.079	22:00	CALIBRATION CHECK RESULT 2		
Air Blank	0.000	22:01	0.079		
Pass			CALIBRATION CHECK RESULT 3		
			0.079		
			MAXIMUM DEVIATION (MUST BE WITHIN 5%)	SPREAD (MUST BE .005 OR LESS)	
			1.2%	0.000	

DIAGNOSTIC TEST RESULTS			RFI TEST RESULTS		
			Test	g/210L	Time
Voltage/Current Test	Pass		Air Blank	0.000	22:01
RAM Test	Pass		Subject Test	RFI*	22:02
EEPROM Checksum Test	Pass		Air Blank	0.000	22:02
Real Time Clock Test	Pass		*RFI Detect		
DSP Test	Pass				
Analytical Stability Test	Pass				
Modem Test	Pass				
Temperature Regulation Test	Pass				
Pass			Pass		

NUMBER OF REFUSALS AND SUBJECT BREATH TESTS IN EACH RANGE SINCE LAST MAINTENANCE REPORT						
REFUSALS	.00-.04	.05-.09	.10-.14	.15-.19	OVER .19	
5	6	3	6	6	11	

List any new parts and describe any alteration or modification that was made to restore the instrument to operate satisfactorily and within established limits (use other side if necessary).

TESTED AND CERTIFIED WITHIN THE DEPARTMENT OF HEALTH AND SENIOR SERVICES GUIDLINES
 TESTED AND CERTIFIED WITHIN THE DEPARTMENT OF HEALTH AND SENIOR SERVICES GUIDLINES

INSPECTING OFFICER		
SIGNATURE <i>Fo Field H. S. Pollard</i>	PRINT NAME POLLARD	
TYPE II PERMIT NUMBER 230167	EXPIRATION DATE 08/14/2015	TELEPHONE NUMBER 816-482-8141



7 Eastgate Dr. • P.O. Box 790 • Jacksonville, IL 62651-0790
217-245-2183 • Fax: 217-243-7634 • www.ilmoproducts.com

Certificate of Analysis

Certificate ID: 6358
Part #: BAC105L080T
Cylinder Size: 105L
Lot Number: 05514080A1
Expiration: 4/1/2016

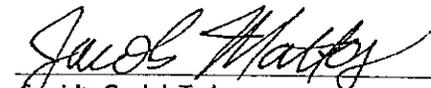
0.080 BAC (For the calibration of instruments used to determine breath alcohol concentration)

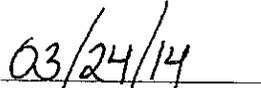
Contents: 105 Liters @ 1000 psig 70°F (21°C)

Component:	Concentration:	Accuracy:	Method:
Ethanol	208 ppm	+/- 0.002 or 2%	NDIR
Nitrogen	Balance	BAC whichever is greater	

*NIST Standard Reference Material
Cylinder No. CCI4290 / Job No. 09160202
Certified 212.8 µmol/mol Ethanol in Nitrogen
for ILMO Products Co., Jacksonville, IL

Store in dry area, away from sources of heat, ignition
and direct sunlight. Do not allow storage area to
exceed 52 °C (125 °F).


Specialty Gas Lab Tech


Date

Distributed by: CMI Inc.
316 East Ninth Street
Owensboro, KY 42303
Phone 866-835-0690
www.alcoholtest.com

