



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES  
 STATE PUBLIC HEALTH LABORATORY  
 BREATH ALCOHOL PROGRAM  
**DATAMASTER MAINTENANCE REPORT**

**RECEIVED**  
 By Carol Day at 10:29 am, Jun 18, 2014

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).  
 Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.  
 Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 990044	NAME OF AGENCY Cape Girardeau Police Department	DATE OF INSPECTION 06/05/2014
LOCATION OF INSTRUMENT (STREET AND CITY) BAT Van		TIME OF INSPECTION 12:50 pm

**CHECKLIST:** Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) 06/05/14 12:50
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER _____ 49 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

INDICATOR LIGHTS

SIMULATOR SOLUTION SUPPLIER RepCo. Marketing Inc. LOT # 13802 EXP. DATE 08/12/2015

SIMULATOR TEMP (34°C ± 0.2°C) 34.0 °C SIMULATOR SN SD2221 EXP. DATE 07/10/2014

CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE  
 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE  
 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 <input checked="" type="checkbox"/> .082	TEST 2 <input checked="" type="checkbox"/> .082	TEST 3 <input checked="" type="checkbox"/> .081
---	---	---

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:  
 (DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	0	(.0-.04)	0	(.05-.09)	1	(.10-.14)	1	(.15-.19)	0	OVER .19	0
----------	---	----------	---	-----------	---	-----------	---	-----------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

**INSPECTING OFFICER**

SIGNATURE 	PRINT FULL NAME Kevin L. Eudy
TYPE II PERMIT NUMBER/EXPIRATION DATE 220132 06/12/2014	TELEPHONE NUMBER (573) 335-6621

RETURN COMPLETED REPORT TO THE:  
 Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office  
 2875 James Blvd.  
 Poplar Bluff, MO 63901

**CERTIFICATE OF ANALYSIS**

**MANUFACTURER AND SUPPLIER: RepCo Marketing, Inc.**  
**LOT NUMBER: 13802**  
**EXPIRATION DATE: August 12, 2015 at 11:59 p.m.**

RepCo Marketing, Inc. certifies the following:

RepCo Marketing, Inc. manufactured, tested and supplied Lot Number 13802 of Alcohol Certified Solution for simulators. Random samples of said lot number were analyzed by an independent laboratory utilizing a gas chromatograph and found to contain .0968 gms/dl +/- .003 gms/dl wt/vol ethanol (95% Confidence).

The alcohol and distilled water used in the solution were found to be free of any interfering substance.

This solution will produce a vapor alcohol value of .080 +/- 3% gms/210L Breath when heated to 34 Degrees Celsius +/- 0.2 Degrees Celsius in a simulator (95% Confidence).

The date of manufacture for this lot number is August 13, 2013.  
The expiration date for this lot number is August 12, 2015 at 11:59 p.m.

This document is a true representation of the original Certificate of Analysis.



Cecil B. Garner, President  
RepCo Marketing, Inc.

STATE OF MASSACHUSETTS  
DEPARTMENT OF HEALTH



PERMIT  
TYPE II



KEVIN LEIDY

is hereby authorized to install and operate equipment, test instruments, inspect, calibrate, repair and service the following brand equipment:

DATABASTER

for the determination of the alcohol content of blood from a sample of expired (intubated) air. Model used for purposes of sections 577.020 through 577.049, RSA 263.

05/12/2012  
220132  
05/12/2014

  
Kevin Leidy  
Director, Department of Public Health

15-000000

15-000000

Face This Side Down - This Edge In First

# BAC DataMaster Evidence Ticket

STATE OF MISSOURI  
SEND LAW ENFORCEMENT ACADEMY

BAC DATAMASTER SERIAL NUMBER 990044  
06/05/14


TESTING OFFICER:  
EUDY/KEVIN/L  
OFFICER I.D.: 186  
PERMIT NUMBER: 220132  
EXPIRATION DATE: 06/12/14  
MISCELLANEOUS DATA:

--- SUPERVISOR MODE ---

BLANK TEST	.000	13:02
INTERNAL STANDARD	VERIFIED	13:03
EXTERNAL STANDARD	.002	13:03
BLANK TEST	.000	13:04
EXTERNAL STANDARD	.002	13:04
BLANK TEST	.000	13:05
EXTERNAL STANDARD	.001	13:05
BLANK TEST	.002	13:06

N = 3  
SIM. = .008  
AVG. = .0016

Operator Signature

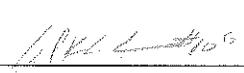

2208-02

Face This Side Down - This Edge In First

# BAC DataMaster Evidence Ticket

STATE OF MISSOURI  
SEND LAW ENFORCEMENT ACADEMY

BAC DATAMASTER SERIAL NUMBER 990044  
06/05/14

ARREST TIME: 12:00  
SUBJECT NAME:  
RFI/TEST  
DOB: 12/23/45 SEX: M  
STATE/D.L.L.: MO/987654  
ARRESTING OFFICER:  
N/R  
OFFICER I.D.: 1  
TESTING OFFICER:  
EUDY/KEVIN/L  
OFFICER I.D.: 186  
PERMIT NUMBER: 220132  
EXPIRATION DATE: 06/12/14  
MISCELLANEOUS DATA:

--- BREATH ANALYSIS ---

BLANK TEST	.000	13:10
INTERNAL STANDARD	VERIFIED	13:10
RADIO INTERFERENCE		

Operator Signature



2208-02

Face This Side Down - This Edge In First

# BAC DataMaster Evidence Ticket

STATE OF MISSOURI  
SEMO LAW ENFORCEMENT ACADEMY

BAC DATAMASTER SERIAL NUMBER 900044  
06/05/14  
12:58

### DIAGNOSTIC CHECK

COMPUTER:	OKAY
PROGRAM (01-97-2000):	OKAY
HEATERS:	
SAMPLE CHAMBER:	OK
FLOW DETECTOR:	OKAY
PUMP	
HIGH SPEED:	OKAY
DETECTOR:	OKAY
FILTERS:	OKAY
QUARTZ STANDARD:	OKAY
CALIBRATION:	OKAY

### PRINTER TEST

! " # \$ % ' ( ) \* + , - . / 0 1 2 3 4 5 6 7 8 9 : ; ( = ) ^ \_ ` a b c d e f g h i j k l m n o p q r s t u v w x y z [ \ ] ^ \_ ` a b c d e f g h i j k l m n o p q r s t u v w x y z ( ) \* +

Operator Signature

