



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
 STATE PUBLIC HEALTH LABORATORY
 BREATH ALCOHOL PROGRAM

RECEIVED
 By Carol Day at 8:56 am, Jan 04, 2016

INTOX EC/IR II MAINTENANCE REPORT

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days). Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service. Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

INTOX EC/IR II SN 12683	NAME OF AGENCY SLMPD	DATE OF INSPECTION 12/31/2015
LOCATION OF INSTRUMENT (STREET AND CITY) 5120 CLAYTON RD ST LOUIS		TIME OF INSPECTION 19:29 CST

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or is operating within established limits. (Write in observed values where determined). Unmarked items must be corrected before using instrument.

DIAGNOSTIC RECORD

<input checked="" type="checkbox"/> BLANK CHECK	<input checked="" type="checkbox"/> CO2 CHECK
<input checked="" type="checkbox"/> FC 1 TEMP	<input checked="" type="checkbox"/> FLOW CHECK
<input checked="" type="checkbox"/> SRC TEMP	<input checked="" type="checkbox"/> FCB CHECK
<input checked="" type="checkbox"/> DET TEMP	<input checked="" type="checkbox"/> CRC COMP CHECK
<input checked="" type="checkbox"/> BT TEMP	<input checked="" type="checkbox"/> CRC CAL CHECK
<input checked="" type="checkbox"/> STD 2 TEMP	<input checked="" type="checkbox"/> PRINT TEST
<input checked="" type="checkbox"/> ETH CHECK	

BREATH ANALYZER ACCURACY STANDARDS

<input type="checkbox"/> SIMULATOR SOLUTION	<input checked="" type="checkbox"/> COMPRESSED ETHANOL-GAS MIXTURE
<input checked="" type="checkbox"/> STANDARD SUPPLIER	INTOXIMETERS
LOT#	AG531302
EXP. DATE	11/09/2017
<input type="checkbox"/> SIMULATOR TEMP (34°C ±0.2°C)	SIMULATOR S/N
	SIMULATOR EXP DATE

CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within +5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

<input type="checkbox"/> 0.10% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
<input checked="" type="checkbox"/> 0.08% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
<input type="checkbox"/> 0.04% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 @ 0.079 g/210L	TEST 2 @ 0.079 g/210L	TEST 3 @ 0.079 g/210L
-----------------------	-----------------------	-----------------------

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:

REFUSALS	6	0-.04	0	.05-.09	0	.10-.14	1	.15-.19	1	OVER .19	2
----------	---	-------	---	---------	---	---------	---	---------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

INSPECTING OFFICER

SIGNATURE 	PRINT FULL NAME CHRISTIAN, SCOTT
TYPE II PERMIT NUMBER 250291	TELEPHONE NUMBER (314) 444-5345
EXPIRATION DATE 12/02/2017	

RETURN COMPLETED REPORT TO THE:
 Breath Alcohol Program, Missouri Department of Health and Senior Services,
 Southeast District Office, 2875 James Blvd, Poplar Bluff, MO 63901



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
STATE PUBLIC HEALTH LABORATORY
BREATH ALCOHOL PROGRAM

RECEIVED

By Carol Day at 8:56 am, Jan 04, 2016

INTOX EC/IR II MAINTENANCE REPORT

REPORT #3

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days). Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service. Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

INTOX EC/IR II SN 12680	NAME OF AGENCY SLMPD BAT VAN	DATE OF INSPECTION 12/31/2015
LOCATION OF INSTRUMENT (STREET AND CITY) 2150 S 59th St ST LOUIS		TIME OF INSPECTION 19:03 CST

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or is operating within established limits. (Write in observed values where determined). Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC RECORD	
<input checked="" type="checkbox"/> BLANK CHECK	<input checked="" type="checkbox"/> CO2 CHECK
<input checked="" type="checkbox"/> FC 1 TEMP	<input checked="" type="checkbox"/> FLOW CHECK
<input checked="" type="checkbox"/> SRC TEMP	<input checked="" type="checkbox"/> PCB CHECK
<input checked="" type="checkbox"/> DET TEMP	<input checked="" type="checkbox"/> CRC COMP CHECK
<input checked="" type="checkbox"/> BT TEMP	<input checked="" type="checkbox"/> CRC CAL CHECK
<input checked="" type="checkbox"/> STD 2 TEMP	<input checked="" type="checkbox"/> PRINT TEST
<input checked="" type="checkbox"/> ETH CHECK	

BREATH ANALYZER ACCURACY STANDARDS

<input type="checkbox"/> SIMULATOR SOLUTION	<input checked="" type="checkbox"/> COMPRESSED ETHANOL-GAS MIXTURE
<input checked="" type="checkbox"/> STANDARD SUPPLIER	INTOXIMETERS
LOT#	AG531302
EXP. DATE	11/09/2017
<input type="checkbox"/> SIMULATOR TEMP (34°C ±0.2°C)	SIMULATOR S/N
	SIMULATOR EXP DATE

CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within +5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

0.10% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE

0.08% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE

0.04% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 ^g 0.079 g/210L	TEST 2 ^g 0.079 g/210L	TEST 3 ^g 0.079 g/210L
----------------------------------	----------------------------------	----------------------------------

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:

REFUSALS	1	0-.04	1	.05-.09	0	.10-.14	0	.15-.19	0	OVER .19	0
----------	---	-------	---	---------	---	---------	---	---------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

INSPECTING OFFICER	
SIGNATURE <i>Scott Christian</i>	PRINT FULL NAME CHRISTIAN, SCOTT
TYPE OF PERMIT NUMBER 250291	EXPIRATION DATE 12/02/2017
	TELEPHONE NUMBER (314) 444-5345

RETURN COMPLETED REPORT TO THE:
Breath Alcohol Program, Missouri Department of Health and Senior Services,
Southeast District Office, 2875 James Blvd, Poplar Bluff, MO 63901



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
 STATE PUBLIC HEALTH LABORATORY
 BREATH ALCOHOL PROGRAM

RECEIVED
 By Carol Day at 8:56 am, Jan 04, 2016

INTOX EC/IR II MAINTENANCE REPORT

REPORT #3

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days). Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service. Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

INTOX EC/IR II SN 12681	NAME OF AGENCY SLMPD	DATE OF INSPECTION 12/31/2015
LOCATION OF INSTRUMENT (STREET AND CITY) 1915 OLIVE ST LOUIS		TIME OF INSPECTION 18:56 CST

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or is operating within established limits. (Write in observed values where determined). Unmarked items must be corrected before using instrument.

DIAGNOSTIC RECORD

<input checked="" type="checkbox"/> BLANK CHECK	<input checked="" type="checkbox"/> CO2 CHECK
<input checked="" type="checkbox"/> FC 1 TEMP	<input checked="" type="checkbox"/> FLOW CHECK
<input checked="" type="checkbox"/> SRC TEMP	<input checked="" type="checkbox"/> FCB CHECK
<input checked="" type="checkbox"/> DET TEMP	<input checked="" type="checkbox"/> CRC COMP CHECK
<input checked="" type="checkbox"/> BT TEMP	<input checked="" type="checkbox"/> CRC CAL CHECK
<input checked="" type="checkbox"/> STD 2 TEMP	<input checked="" type="checkbox"/> PRINT TEST
<input checked="" type="checkbox"/> ETH CHECK	

BREATH ANALYZER ACCURACY STANDARDS

<input type="checkbox"/> SIMULATOR SOLUTION	<input checked="" type="checkbox"/> COMPRESSED ETHANOL-GAS MIXTURE		
<input checked="" type="checkbox"/> STANDARD SUPPLIER	INTOXIMETERS	LOT# AG531302	EXP. DATE 11/09/2017
<input type="checkbox"/> SIMULATOR TEMP (34°C ±0.2°C)	SIMULATOR S/N	SIMULATOR EXP DATE	

CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

<input type="checkbox"/> 0.10% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
<input checked="" type="checkbox"/> 0.08% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
<input type="checkbox"/> 0.04% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 0.080 g/210L	TEST 2 0.080 g/210L	TEST 3 0.080 g/210L
---------------------	---------------------	---------------------

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:

REFUSALS	3	0-.04	0	.05-.09	0	.10-.14	0	.15-.19	0	OVER .19	0
----------	---	-------	---	---------	---	---------	---	---------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

INSPECTING OFFICER

SIGNATURE <i>Scott Christian</i>	PRINT FULL NAME CHRISTIAN, SCOTT
TYPE II PERMIT NUMBER 250291	EXPIRATION DATE 12/02/2017
	TELEPHONE NUMBER (314) 444-5345

RETURN COMPLETED REPORT TO THE:
 Breath Alcohol Program, Missouri Department of Health and Senior Services,
 Southeast District Office, 2875 James Blvd, Poplar Bluff, MO 63901



Airgas USA LLC (LAB)
 3500 Bernard Street
 St. Louis, Mo. 63103
 Ph: (314) 533-3100
 Fax: (314) 533-7328

Certificate of Analysis

Customer Name
 Exclusive Supplier
 Intoximeters, Inc.
 2081 Craig Road
 St. Louis, Mo 63146

Test Date: 10-Nov-2015

Lot # AG531302 Model 108cacd


<u>Exp. Date</u>	<u>Cyl. Type</u>	<u>Component</u>	<u>Certified Concentration</u>
9-Nov-2017	108	Ethanol Nitrogen	0.080 ± 0.002 BrAC (218 ppm) Balance

Certification Traceable to N.I.S.T. RGM Ethanol Standards:

<u>Serial No.</u>	<u>Concentration</u>	<u>Serial No.</u>	<u>Concentration</u>
EB0010581	391.8 ppm	EB0010603	392.5 ppm
EB0010570	259.8 ppm	EB0010559	258.9 ppm
EB0010285	209.0 ppm	EB0010595	208.9 ppm
EB0010561	103.7 ppm	EB0010562	104.9 ppm
EB0010681	52.22 ppm	EB0010579	52.94 ppm

Analytical Method: NDIR


Digitally signed by Quality Control
 Date: 2015.11.10 09:43:23 -06:00
 Reason: Dry gas standard certification of analysis
 Location: Airgas USA LLC (Lab)

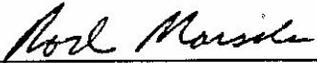Analyst: 
 Rod Marsala

ISO 17025:2005 A2LA accredited. Certificate Number 2989.01



STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

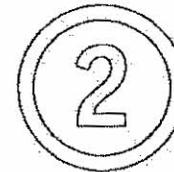

**PERMIT
 TYPE II**

SCOTT CHRISTIAN

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

INTOX EC/IR II

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 12/2/2015

NUMBER 250291

EXPIRES 12/2/2017

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES