

RECEIVED

By Carol Day at 3:23 pm, Feb 19, 2016



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
STATE PUBLIC HEALTH LABORATORY
BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

REPORT #6

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).
Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.
Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 201203	NAME OF AGENCY Chaffee Police Department	DATE OF INSPECTION 06-22-2015
LOCATION OF INSTRUMENT (STREET AND CITY) 222 W. Yorkum Ave Chaffee		TIME OF INSPECTION 08:19 hours

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) 06-22-2015 @ 08:19 hours
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER +49 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

<input checked="" type="checkbox"/> INDICATOR LIGHTS
<input checked="" type="checkbox"/> SIMULATOR SOLUTION SUPPLIER RepCo Marketing LOT # 14001 EXP. DATE 4-30-2016
<input checked="" type="checkbox"/> SIMULATOR TEMP (34°C ± 0.2°C) +34.0 °C SIMULATOR SN SD2228 EXP. DATE 3-12-2016

CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 • .097	TEST 2 • .097	TEST 3 • .097
----------------------	----------------------	----------------------

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:
(DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS 0	(0-.04) 1	(.05-.09) 0	(.10-.14) 0	(.15-.19) 0	OVER .19 0
-------------------	------------------	--------------------	--------------------	--------------------	-------------------

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).
None.

INSPECTING OFFICER

SIGNATURE William Sammut	PRINT FULL NAME William Sammut
TYPE II PERMIT NUMBER/EXPIRATION DATE 250014 1-14-2017	TELEPHONE NUMBER (573) 887-6911

RETURN COMPLETED REPORT TO THE: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
2875 James Blvd.
Poplar Bluff, MO 63901

Face This Side Down - This Edge In First

BAC DataMaster
Evidence Ticket

STATE OF MISSOURI
CHAFFEE POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 201203
06/22/15
08:19

--- DIAGNOSTIC CHECK ---

COMPUTER: OKAY
PROGRAM (04-07-2009): OKAY
HEATERS
SAMPLE CHAMBER: 49c
FLOW DETECTOR: OKAY
PUMP
HIGH SPEED: OKAY
DETECTOR: OKAY
FILTERS: OKAY
QUARTZ STANDARD: OKAY
CALIBRATION: OKAY

PRINTER TEST

!"#\$%&'()*+,-./0123456789:;<=>?@ABCDEFGHIJ
KLMNOPQRSTUVWXYZ[\]^_`abcdefghijklmnopqrstuvwxyz{|}~
!@#\$%^&*()~

Operator Signature

W. Summit SEMO
OPS 521

Face This Side Down - This Edge In First

BAC DataMaster
Evidence Ticket

STATE OF MISSOURI
CHAFFEE POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 201203
06/22/15

TESTING OFFICER:
SAMMUT/WILLIAM/J
OFFICER I.D.: SEMO 521
PERMIT NUMBER: 250014
EXPIRATION DATE: 01/14/17
MISCELLANEOUS DATA:
JUNE 2015 MAINT REPORT

--- SUPERVISOR MODE ---

BLANK TEST	.000	09:13
INTERNAL STANDARD	VERIFIED	09:13
EXTERNAL STANDARD	.097	09:13
BLANK TEST	.000	09:14
EXTERNAL STANDARD	.097	09:14
BLANK TEST	.000	09:15
EXTERNAL STANDARD	.097	09:15
BLANK TEST	.000	09:16

N = 3
SIM. = .1
AVG. = .097

Operator Signature

W. Sammut SEMO
DPS 521



STATE OF MISSOURI
DEPARTMENT OF HEALTH AND SENIOR SERVICES
BREATH ALCOHOL PROGRAM



PERMIT
TYPE II

WILLIAM J SAMMUT

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

DATAMASTER

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 1/14/2015
NUMBER 250014
EXPIRES 1/14/2017

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

Darl Vesterly
acting director

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

MO-580-0771 (6-10)

LAB-4 (R6-10)

CERTIFICATE OF ANALYSIS

MANUFACTURER AND SUPPLIER: RepCo Marketing, Inc.

LOT NUMBER: 14001

EXPIRATION DATE: April 30, 2016 at 11:59 p.m.

RepCo Marketing, Inc. certifies the following:


RepCo Marketing, Inc. manufactured, tested and supplied Lot Number 14001 of Alcohol Certified Solution for simulators. Random samples of said lot number were analyzed by an independent laboratory utilizing a gas chromatograph and found to contain .1216 gms/dl +/- .003 gms/dl wt/vol ethanol (95% Confidence).

The alcohol and distilled water used in the solution were found to be free of any interfering substance.

This solution will produce a vapor alcohol value of .100 +/- 3% gms/210L Breath when heated to 34 Degrees Celsius +/- 0.2 Degrees Celsius in a simulator (95% Confidence).

The date of manufacture for this lot number is May 1, 2014
The expiration date for this lot number is April 30, 2016 at
11:59 p.m.

This document is a true representation of the original Certificate of Analysis.

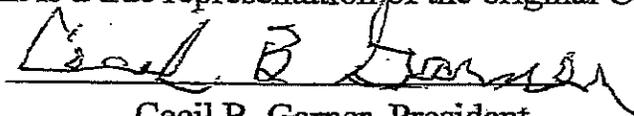

Cecil B. Garner, President
RepCo Marketing, Inc.