



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
 STATE PUBLIC HEALTH LABORATORY
 BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

RECEIVED REPORT #6
 By Carol Day at 8:57 am, Mar 10, 2015

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).
 Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.
 Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 204112	NAME OF AGENCY MSHP	DATE OF INSPECTION 03/03/2015
LOCATION OF INSTRUMENT (STREET AND CITY) RANDOLPH CO SO 372 HWY JJ HUNTSVILLE		TIME OF INSPECTION 8:50 pm

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) 03/03/2015 2059
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 50°C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

<input checked="" type="checkbox"/> INDICATOR LIGHTS
<input checked="" type="checkbox"/> SIMULATOR SOLUTION SUPPLIER Guth LOT # 13290 EXP. DATE 10/29/2015
<input checked="" type="checkbox"/> SIMULATOR TEMP (34°C ± 0.2°C) 34.0 °C SIMULATOR SN MP2518 EXP. DATE 09/19/2015

CALIBRATION CHECK -- (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 • 100%	TEST 2 • 101%	TEST 3 • 102%
---------------	---------------	---------------

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

**INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:
 (DO NOT INCLUDE SELF-ADMINISTERED TESTS)**

REFUSALS	0	(0-.04)	0	(.05-.09)	0	(.10-.14)	2	(.15-.19)	1	OVER .19	1
----------	---	---------	---	-----------	---	-----------	---	-----------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

INSTRUMENT CALIBRATED

INSPECTING OFFICER	
SIGNATURE 	PRINT FULL NAME Bruce E. Thoman
TYPE II PERMIT NUMBER/EXPIRATION DATE 240185 04/22/2016	TELEPHONE NUMBER (660) 385-2132

RETURN COMPLETED REPORT TO THE:
 Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
 2875 James Blvd.
 Poplar Bluff, MO 63901



GUTH LABORATORIES, INC.

889 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-664-5470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 13290 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on October 31, 2013, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1202% (w/vol) ethyl alcohol. The expiration date for this lot number is October 29, 2015 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

**Ted L. Pauley, President
GUTH LABORATORIES, INC.**

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN122211-02 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

BAC DataMaster
Evidence Ticket

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204112
03/03/15

ARREST TIME: 18:00
SUBJECT NAME: RFI/TEST
DOB: 11/11/11 SEX: M
STATE/D.L.: MO/111111
ARRESTING OFFICER:
X
OFFICER I.D.: X
TESTING OFFICER:
THOMAN/BRUCE/E
OFFICER I.D.: 1256
PERMIT NUMBER: 240185
EXPIRATION DATE: 04/22/16
MISCELLANEOUS DATA:

--- BREATH ANALYSIS ---

BLANK TEST .000 21:11
INTERNAL STANDARD VERIFIED 21:11
RADIO INTERFERENCE

Operator Signature



BAC DataMaster
Evidence Ticket


MISSOURI STATE HIGHWAY PATROL
BAC DATA MASTER SERIAL NUMBER 204112
03/03/15

TESTING OFFICER:
THOMAS/BRUCE/E
OFFICER I.D.: 1256
PERMIT NUMBER: 240185
EXPIRATION DATE: 04/22/16
MISCELLANEOUS DATA:

--- SUPERVISOR MODE ---

BLANK TEST	.000	21:04
INTERNAL STANDARD	VERIFIED	21:05
EXTERNAL STANDARD	.100	21:05
BLANK TEST	.000	21:05
EXTERNAL STANDARD	.101	21:06
BLANK TEST	.000	21:07
EXTERNAL STANDARD	.102	21:07
BLANK TEST	.000	21:07

N = 3
SIM. = .101
AVG. = .101

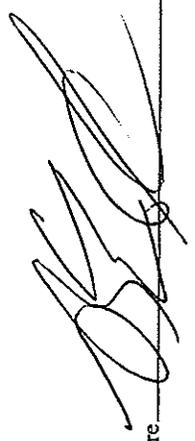

Operator Signature

BAC DataMaster
Evidence Ticket

MISSOURI STATE HIGHWAY PATROL
BAC DATA MASTER SERIAL NUMBER 204112
03/03/15
20:59

--- DIAGNOSTIC CHECK ---

COMPUTER: OKAY

PROGRAM (04-07-2009): OKAY

HEATERS

SAMPLE CHAMBER: 50C

FLOW DETECTOR: OKAY

PUMP

HIGH SPEED: OKAY

DETECTOR: OKAY


FILTERS: OKAY

QUARTZ STANDARD: OKAY

CALIBRATION: OKAY

PRINTER TEST

!"#\$%&'()*+,-./0123456789:;<=>?@ABCDEFGHIJKLmnop
HIJKLMNOPQRSTUVWXYZ[\]^_`abcdefgijklmnop
qrs'tuvwxyz{|}~

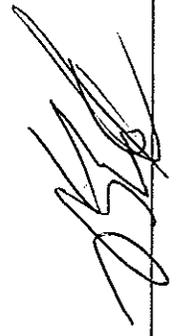

Operator Signature



STATE OF MISSOURI
DEPARTMENT OF HEALTH AND SENIOR SERVICES
BREATH ALCOHOL PROGRAM

2

**PERMIT
TYPE II**

BRUCE E THOMAN

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

DATAMASTER, INTOX DMT

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 4/22/2014

NUMBER 240185

EXPIRES 4/22/2016

MO 590-0771 (6-10)

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

LAB-4 (R6-10)

 STATE OF MISSOURI
DEPARTMENT OF HEALTH AND SENIOR SERVICES
BREATH ALCOHOL PROGRAM

INSTRUMENT OPERATOR CARD

The named cardholder is authorized to operate an evidential breath alcohol instrument for the determination of the alcoholic content in breath form of expired air in Missouri.



Operator THOMAN, BRUCE
Permit No 240185
Date Issued 4/22/2014 Date Expires 4/22/2016