



BREATH-ALCOHOL PROGRAM  
**DATAMASTER MAINTENANCE REPORT**

REPORT #6

**RECEIVED**

By Carol Day at 9:14 am, Nov 17, 2015

Complete this report at the time of the regular monthly preventive maintenance of the instrument. Complete this report whenever the instrument is serviced or repaired and whenever the instrument is used. Retain the original and send a copy within 15 days to the Breath Alcohol Program.

DATAMASTER SN 960062	NAME OF AGENCY Riverview PD	DATE OF INSPECTION 11/11/2015
LOCATION OF INSTRUMENT (STREET AND CITY) 9699 Lilac Riverview		TIME OF INSPECTION 9:45 am

**CHECKLIST:** Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) 11/11/2015 10:04
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 50 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

INDICATOR LIGHTS

SIMULATOR SOLUTION SUPPLIER Guth LOT # 15120 EXP. DATE 04/29/2017

SIMULATOR TEMP (34°C ± 0.2°C) 34 °C SIMULATOR SN SD2769 EXP. DATE 06/26/2016

CALIBRATION CHECK – (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE  
 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE  
 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 <u>.099</u>	TEST 2 <u>.099</u>	TEST 3 <u>.099</u>
--------------------	--------------------	--------------------

PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT: (DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	0	(0-.04)	0	(.05-.09)	0	(.10-.14)	0	(.15-.19)	0	OVER .19	0
----------	---	---------	---	-----------	---	-----------	---	-----------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

<b>INSPECTING OFFICER</b>	
SIGNATURE 	PRINT FULL NAME Nicholes C. Allred
TYPE II PERMIT NUMBER/EXPIRATION DATE 250001 01/02/2017	TELEPHONE NUMBER (314) 868-9130

RETURN COMPLETED REPORT TO THE: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office  
 2875 James Blvd.  
 Poplar Bluff, MO 63901



**GUTH LABORATORIES, INC.**

690 NORTH 67th STREET • HARRISBURG, PA 17111-4611 • TELEPHONE: 717-564-5470

**CERTIFICATE OF ANALYSIS**

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 15120 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on May 4, 2015, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1209% (w/vol) ethyl alcohol. The expiration date for this lot number is April 29, 2017 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President  
GUTH LABORATORIES, INC.

***NIST Traceability:***

*Testing was conducted using Cerilliant Reference Standard lot number FN08051301 whose values are traceable to NIST.*

*All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.*

BAC DataMaster  
Evidence Ticket

STATE OF MISSOURI  
RIVERVIEW POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 360062  
11/11/15  
13:04

--- DIAGNOSTIC CHECK ---

COMPUTER: OKAY  
PROGRAM (84-97 2000): OKAY  
SERTERS  
SAMPLE CHAMBER: 300  
LOW DETECTOR: OKAY  
MMP  
HIGH SPEED: OKAY  
DETECTOR: OKAY  
FILTERS: OKAY  
QUARTZ STANDARD: OKAY  
CALIBRATION: OKAY

PRINTER TEST

!"#\$%&'()\*+,-./0123456789:;<=>?@ABCDEFGHIJKL  
MNOPQRSTUVWXYZ[\]^\_`abc`defghijklmnop  
qrstuvwxyz{|}~\*+</p>
</div>
<div data-bbox="15 729 87 742" data-label="Text">

Signature

2208-02

BAC DataMaster  
Evidence Ticket

STATE OF MISSOURI  
RIVERVIEW POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 360062  
11/11/15

ARREST TIME: 09:45  
SUBJECT NAME:  
TEST: JOP/JAY  
DOB: 11/22/66 SEX: M  
STATE/D.L.# MO/YT/PTRE  
ARRESTING OFFICER:  
ALLRED/NICHOLS/C  
OFFICER I.D.: 247  
TESTING OFFICER:  
ALLRED/NICHOLS/C  
OFFICER I.D.: 247  
PERMIT NUMBER: 250001  
EXPIRATION DATE: 01/02/17  
MISCELLANEOUS DATA:

--- BREATH ANALYSIS ---

BLANK TEST .000 18:16  
INTERNAL STANDARD VERIFIED 19:16  
SUBJECT SAMPLE .000 18:17  
RADIO INTERFERENCE

Operator Signature

2209-02

# BAC DataMaster Evidence Ticket

STATE OF MISSOURI  
RIVERVIEW POLICE DEPARTMENT

BAC DATA MASTER SERIAL NUMBER 960763  
11/11/15


TESTING OFFICER:  
ALLRED/MICHAEL/C  
OFFICER I.D. # 247  
PERMIT NUMBER: 250001  
EXPIRATION DATE: 01-03-17  
MISCELLANEOUS DATA:

--- SUPERVISOR CODE ---

BLANK TEST	.000	13:00
INTERNAL STANDARD	VERIFIED	13:02
EXTERNAL STANDARD	.000	13:03
BLANK TEST	.000	13:10
EXTERNAL STANDARD	.000	13:13
BLANK TEST	.000	13:11
EXTERNAL STANDARD	.000	13:11
BLANK TEST	.000	13:12

N = 3  
STDEV. = .1  
AVE. = .000

Operator Signature

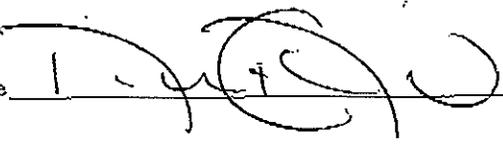



STATE OF MISSOURI  
 DEPARTMENT OF HEALTH AND SENIOR SERVICES  
 BREATH ALCOHOL PROGRAM



**PERMIT**  
**TYPE II**  
**NICHOLE S C ALLRED**

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

**DATAMASTER**

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 1/2/2015

NUMBER 250001

EXPIRES 1/2/2017

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

acting director

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

MO 580-0771 (6-10)

LAB-4 (R9-10)

STATE OF MISSOURI  
 DEPARTMENT OF HEALTH AND SENIOR SERVICES  
 BREATH ALCOHOL PROGRAM

**INSTRUMENT OPERATOR CARD**

*The named cardholder is authorized to operate an evidential breath alcohol instrument for the determination of the alcoholic content in breath form of expired air in Missouri.*

Operator ALLRED, NICHOLE S  
 Permit No 250001  
 Date Issued 1/2/2015 Date Expires 1/2/2017