

RECEIVED

By Stephen Wilson at 10:01 am, May 27, 2020



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
STATE PUBLIC HEALTH LABORATORY
BREATH ALCOHOL PROGRAM

INTOX EC/IR II MAINTENANCE REPORT

REPORT #3

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days). Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service. Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

INTOX EC/IR II SN 12679	NAME OF AGENCY Chesterfield PD	DATE OF INSPECTION 05/26/2020
----------------------------	-----------------------------------	----------------------------------

LOCATION OF INSTRUMENT (STREET AND CITY) 690 Chesterfield Pkwy W Chesterfield	TIME OF INSPECTION 10:21 CST
--	---------------------------------

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or is operating within established limits. (Write in observed values where determined). Unmarked items must be corrected before using instrument.

 DIAGNOSTIC RECORD

<input checked="" type="checkbox"/> BLANK CHECK	<input checked="" type="checkbox"/> CO2 CHECK
<input checked="" type="checkbox"/> FC 1 TEMP	<input checked="" type="checkbox"/> FLOW CHECK
<input checked="" type="checkbox"/> SRC TEMP	<input checked="" type="checkbox"/> FCB CHECK
<input checked="" type="checkbox"/> DET TEMP	<input checked="" type="checkbox"/> CRC COMP CHECK
<input checked="" type="checkbox"/> BT TEMP	<input checked="" type="checkbox"/> CRC CAL CHECK
<input checked="" type="checkbox"/> STD 2 TEMP	<input checked="" type="checkbox"/> PRINT TEST
<input checked="" type="checkbox"/> ETH CHECK	

BREATH ANALYZER ACCURACY STANDARDS

<input type="checkbox"/> SIMULATOR SOLUTION	<input checked="" type="checkbox"/> COMPRESSED ETHANOL-GAS MIXTURE
<input checked="" type="checkbox"/> STANDARD SUPPLIER ILMO	LOT# 1781908011 EXP. DATE 08/05/2021
<input type="checkbox"/> SIMULATOR TEMP (34°C ±0.2°C)	SIMULATOR S/N SIMULATOR EXP DATE

 CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within $\pm 5\%$ of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

<input type="checkbox"/> 0.10% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
<input checked="" type="checkbox"/> 0.08% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
<input type="checkbox"/> 0.04% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 → 0.079 g/210L	TEST 2 → 0.079 g/210L	TEST 3 → 0.079 g/210L
-----------------------	-----------------------	-----------------------

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:

REFUSALS	0	0-.04	19	.05-.09	0	.10-.14	2	.15-.19	0	OVER .19	0
----------	---	-------	----	---------	---	---------	---	---------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

INSTRUMENT IS OPERATING WITHIN MO DEPT OF HEALTH SPECIFICATIONS

Instrument
G.R. #264

INSPECTING OFFICER

SIGNATURE ▶ <i>G. Rupp #264</i>	PRINT FULL NAME RUPP, GREGORY
TYPE II PERMIT NUMBER 290094	TELEPHONE NUMBER (636) 537-3000
EXPIRATION DATE 04/22/2021	

RETURN COMPLETED REPORT TO THE:

Breath Alcohol Program, Missouri Department of Health and Senior Services,
Southeast District Office, 2875 James Blvd, Poplar Bluff, MO 63901



7 Eastgate Dr. • P.O. Box 790 • Jacksonville, IL 62651-0790
217-245-2183 • Fax: 217-243-7634 • www.ilmoproducts.com

Certificate of Analysis

Certificate ID: 12199
Part #: BAC108L080I
Cylinder Size: 108L
Lot Number: 17819080I1
Expiration: 8/5/2021

0.080 BAC (For the calibration of instruments used to determine breath alcohol concentration)


Contents: 108 Liters @ 1200 psig 70°F (21°C)

Component:	Reported Concentration:	Analytical Accuracy (U, k=2):	Analytical Method:
Ethanol	218 ppm	+/- 0.002 BAC (G/210L) [5.2 ppm]	NDIR
Nitrogen	Balance		


***Traceable to:**

Certified Reference Material - 262.4 µmol/mol
Ethanol in Nitrogen - Serial No. GN0015026 Lot No. 050319E11

Store in dry area, away from sources of heat, ignition and direct sunlight. Do not allow storage area to exceed 52 °C (125 °F).


Specialty Gas Lab Tech

07-05-19
Date

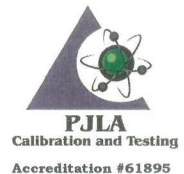

The calibration results within this certificate were obtained using equipment and standards capable of producing analytical results traceable to NIST, and apply only to the items contained on this certificate. ILMO Products Company makes no warranty or representation as to the suitability of the use of any information provided for any particular purpose. The information use is at the sole discretion and risk of the user. Liability shall be limited to established replacement cost of this material or service.

ISO/IEC 17025:2005 Accredited Laboratory



STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM



PERMIT
TYPE II
GREGORY RUPP

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

ALCO-SENSOR IV WITH PRINTER, INTOX EC/IR II

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 4/22/2019

NUMBER 290094

EXPIRES 4/22/2021

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

 **STATE OF MISSOURI**
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

INSTRUMENT OPERATOR CARD

The named cardholder is authorized to operate an evidential breath alcohol instrument for the determination of the alcoholic content in breath form of expired air in Missouri.

Operator RUPP, GREGORY
Permit No 290094
Date Issued 4/22/2019 **Date Expires** 4/22/2021

