



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
 STATE PUBLIC HEALTH LABORATORY
ALCO-SENSOR IV WITH PRINTER MAINTENANCE REPORT

REPORT #7

Complete this report in duplicate at the time of the regular monthly preventative maintenance check, and whenever instrument is repaired. Send copy to Department of Health and Senior Services; retain original in department file.

ALCO SENSOR IV SN 104658	PRINTER SN 096.3580.740	DATE OF INSPECTION 09/03/2022
LOCATION OF INSTRUMENT (STREET AND CITY) 3101 Miller Street, Bethany, MO, 64424		TIME OF INSPECTION 11:10 am

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIGITAL READOUT (ALL ELEMENTS OPERATIONAL)
<input checked="" type="checkbox"/> TEMPERATURE OF ALCO SENSOR (10°C - 40°C)
<input checked="" type="checkbox"/> PRINTER WORKING PROPERLY
<input checked="" type="checkbox"/> TIME AND DATE DISPLAYING PROPERLY

BREATH ALCOHOL ACCURACY STANDARDS

<input checked="" type="checkbox"/> SIMULATOR SOLUTION	<input type="checkbox"/> COMPRESSED ETHANOL-GAS MIXTURE
<input checked="" type="checkbox"/> STANDARD SUPPLIER <u>Repco</u>	LOT # <u>20001</u> EXP. DATE <u>10/07/2022</u>
<input checked="" type="checkbox"/> SIMULATOR TEMPERATURE (34°C ± 0.2°C) <u>34.0</u>	SIMULATOR SN <u>MP2502</u> SIMULATOR EXP DATE <u>02/08/2023</u>
<input checked="" type="checkbox"/> CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT) Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Check the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)	
<input checked="" type="checkbox"/> 0.100% STANDARD - MUST READ BETWEEN 0.095% and 0.105% INCLUSIVE	
<input type="checkbox"/> 0.080% STANDARD - MUST READ BETWEEN 0.076% and 0.084% INCLUSIVE	
<input type="checkbox"/> 0.040% STANDARD - MUST READ BETWEEN 0.038% and 0.042% INCLUSIVE	

TEST 1 ← .099	TEST 2 ← .099	TEST 3 ← .099
---------------	---------------	---------------

RFI DETECTOR OPERATING

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:
 (DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	0	(.0-.04)	0	(.05-.09)	0	(.10-.14)	0	(.15-.19)	0	(OVER .19)	0
----------	---	----------	---	-----------	---	-----------	---	-----------	---	------------	---

List any new parts and describe any alteration or modification that was made to restore the instrument to operate satisfactorily and within established limits (use other side if necessary).

INSPECTING OFFICER	
SIGNATURE 	PRINT NAME Corporal M. J. Miller #1418
TYPE II PERMIT NUMBER/EXPIRATION DATE 200292	TELEPHONE NUMBER (816) 387-2345

Return completed report to the: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
 2875 James Boulevard
 Poplar Bluff, MO 63901



RepCo Marketing Co
3101-188 Stony Brook Drive
Raleigh, NC 27604
919-876-5480

CERTIFICATE OF ANALYSIS

MANUFACTURER AND SUPPLIER: RepCo Marketing Co.
LOT NUMBER: 20001
EXPIRATION DATE: October 7, 2022 at 11:59 p.m.

RepCo Marketing Co. certifies the following:

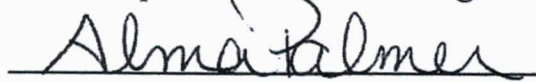
RepCo Marketing Co. prepared, tested and supplied Lot Number 20001 of Alcohol Certified Solution for simulators. This solution was manufactured and tested by RepCo Marketing Co., with confirmation by ISO 17025 and ISO 17034 accredited institution Alcohol Countermeasure Systems, using NIST standards. Random samples were analyzed by Alcohol Countermeasure Systems utilizing a gas chromatograph and found to contain .1227 gms/dl +/- .003 gms/dl wt/vol ethanol (95% Confidence).

The alcohol and distilled water used in the solution were found to be free of any interfering substance.

This solution will produce a vapor alcohol value of .100 +/- 3% gms/210L breath when heated to 34 Degrees Celsius +/- 0.2 Degrees Celsius in a simulator (95% Confidence).

The date of manufacture for this lot number is October 8, 2020 The expiration date for this lot number is October 7, 2022 at 11:59 p.m.

This document is a true representation of the original Certificate of Analysis.


Alma Palmer, Operations Manager
RepCo Marketing Co.



STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

2

PERMIT
TYPE II
MICHAEL J. MILLER

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

ALCO-SENSOR IV WITH PRINTER, INTOX DMT

for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

DATE 12/8/2020

NUMBER 200292

EXPIRES 12/8/2022

MO 580 0771 (6-10)

DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

LAB-4 (R6-10)

STATE OF MISSOURI
 DEPARTMENT OF HEALTH AND SENIOR SERVICES
 BREATH ALCOHOL PROGRAM

INSTRUMENT OPERATOR CARD

The named cardholder is authorized to operate an evidential breath alcohol instrument for the determination of the alcoholic content in breath form of expired air in Missouri

Operator MILLER, MICHAEL
Permit No 200292
Date Issued 12/8/2020 **Date Expires** 12/8/2022